


Alcohols

CUSTOM
DEDICATED

MEGA®
gas chromatography solutions

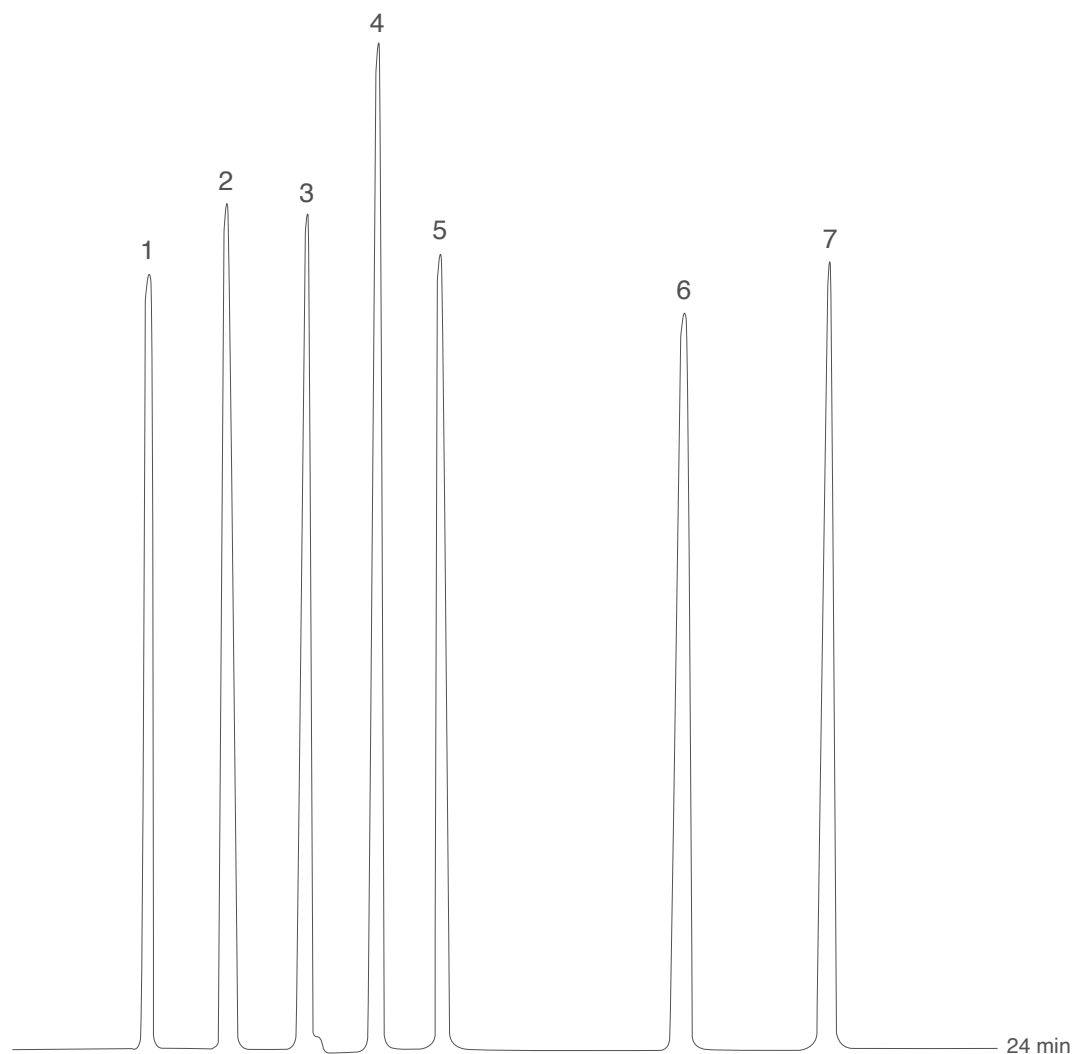

Column: **MEGA-PS255 - 0.32mm, 3.00 μ m, 25m**

Catalog Code: C-PS255-032-300-25

Conditions:

Injection: Split 250°C, 1:50 Split Ratio.

Oven Temp.: 50°C, 3°C/min, 120°C.

Detector: FID, 250°C.

Carrier Gas: Nitrogen, 0.4 mL/min.

Peaks:

1. Methanol
2. Ethanol
3. 2-propanol
4. 1-propanol
5. Tert-butanol
6. Iso-butanol
7. 1-butanol