

(+/-)-2,3-diphenylpropionic acid

CHIRALDEX

FAST

MEGA[®]
gas chromatography solutions

Chiral separation on MEGA-DEX columns


Column: **MEGA-DEX DAC Beta FAST - 0.10mm, 0.10 μ m, 5m**

Catalog Code: DEXF-DAC-B-010-010-5

Conditions:

Oven program: 150°C, isothermal.

Carrier Gas: Helium, 0.4mL/min.

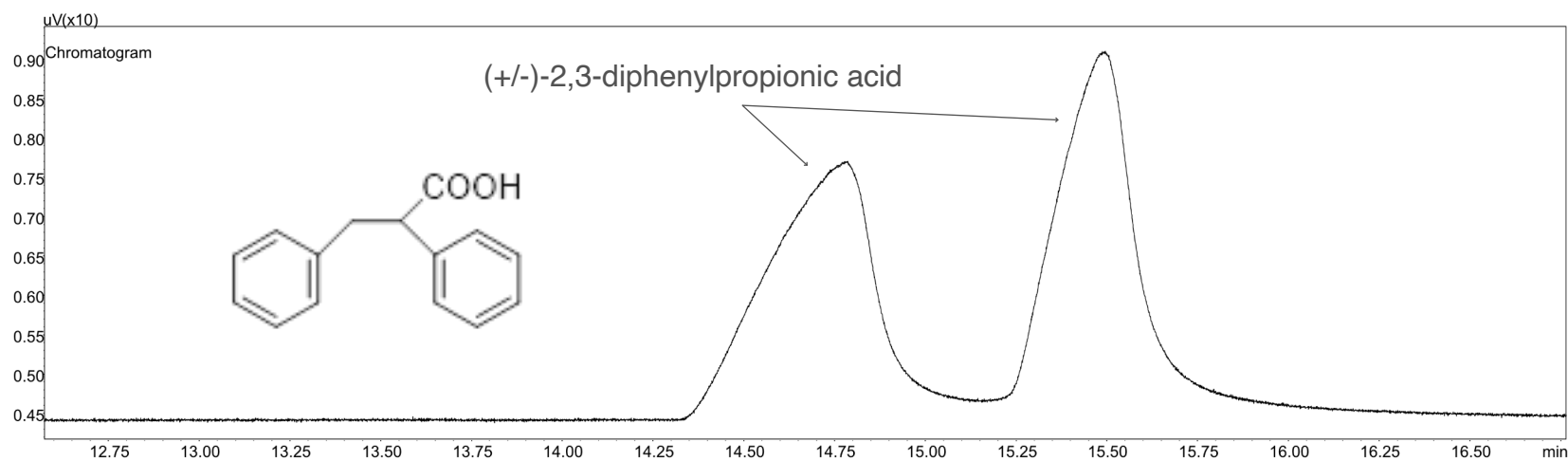

Acknowledgments: many thanks to Barbara Sgorbini, Cecilia Cagliero, Chiara Cordero and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.