

## 2-methyl butyric acid (+/-)


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 $\mu$ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

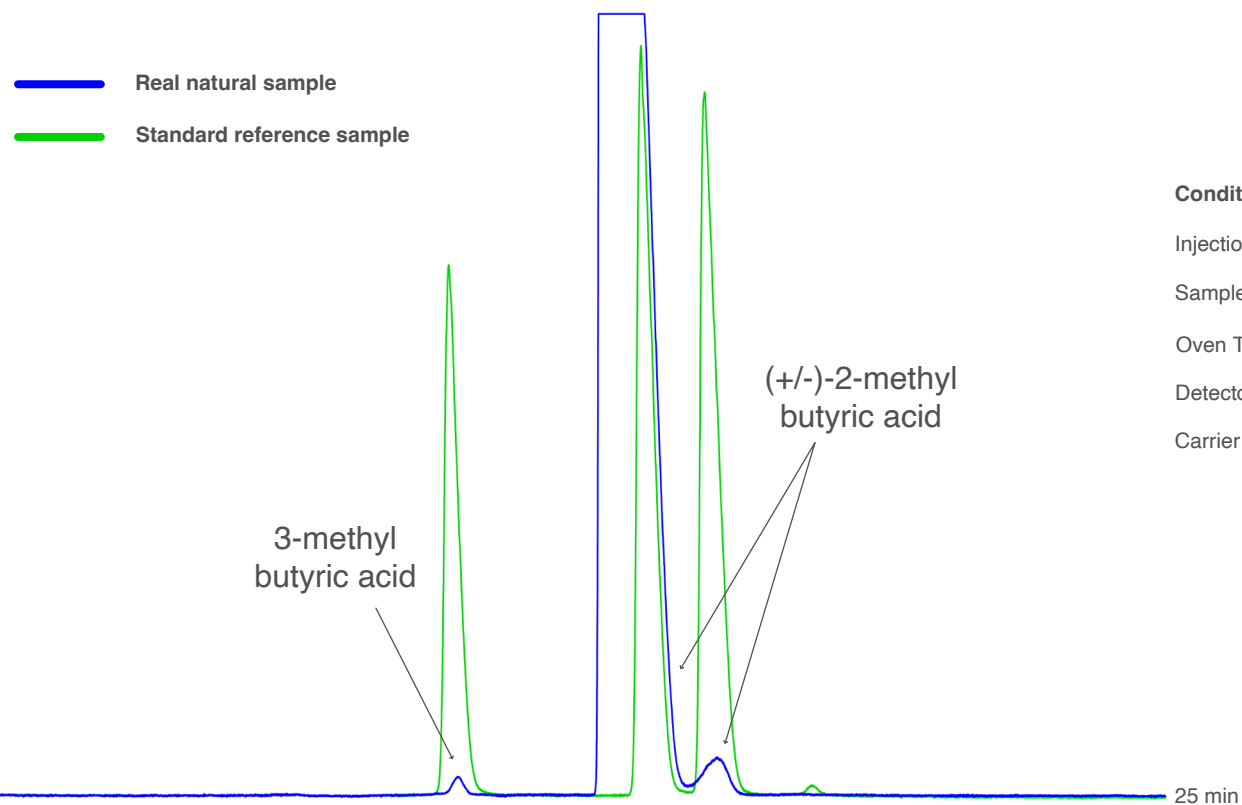

### Conditions:

Injection: Split 250°C, 0.5 $\mu$ L, 350mL/min Split Flow.

Sample dilution: 10 $\mu$ L/1mL each in Methanol.

Oven Temp.: 50°C, 2°C/min.

Detector: FID, 250°C.

Carrier Gas: Hydrogen, 80kPa.