

3-hydroxytetrahydrofuran - enantiomers


CHIRAL DEX

MEGA[®]
gas chromatography solutions

Chiral separation on MEGA-DEX columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

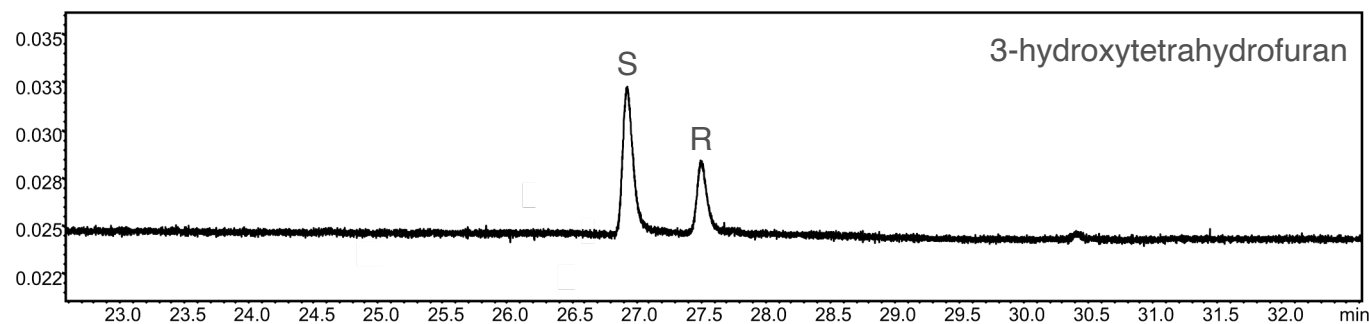


Conditions:

Oven Temp.: 50°C, 2°C/min.

Carrier Gas: Hydrogen, 1mL/min.

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DAC-B-025-025-25

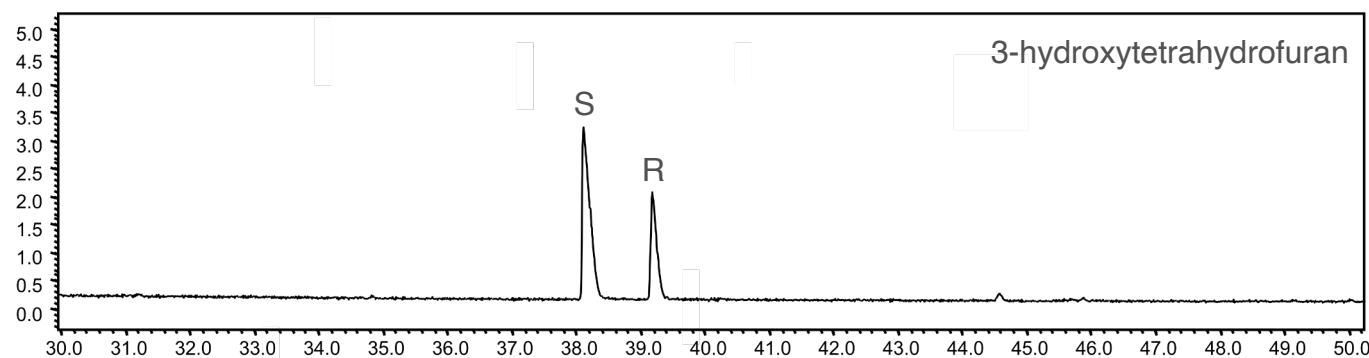

Acknowledgments: many thanks to Barbara Sgorbini, Cecilia Cagliero, Chiara Cordero and Carlo Biochi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.