


Basic Drugs - Underivatized

STANDARD

MEGA[®]
gas chromatography solutions

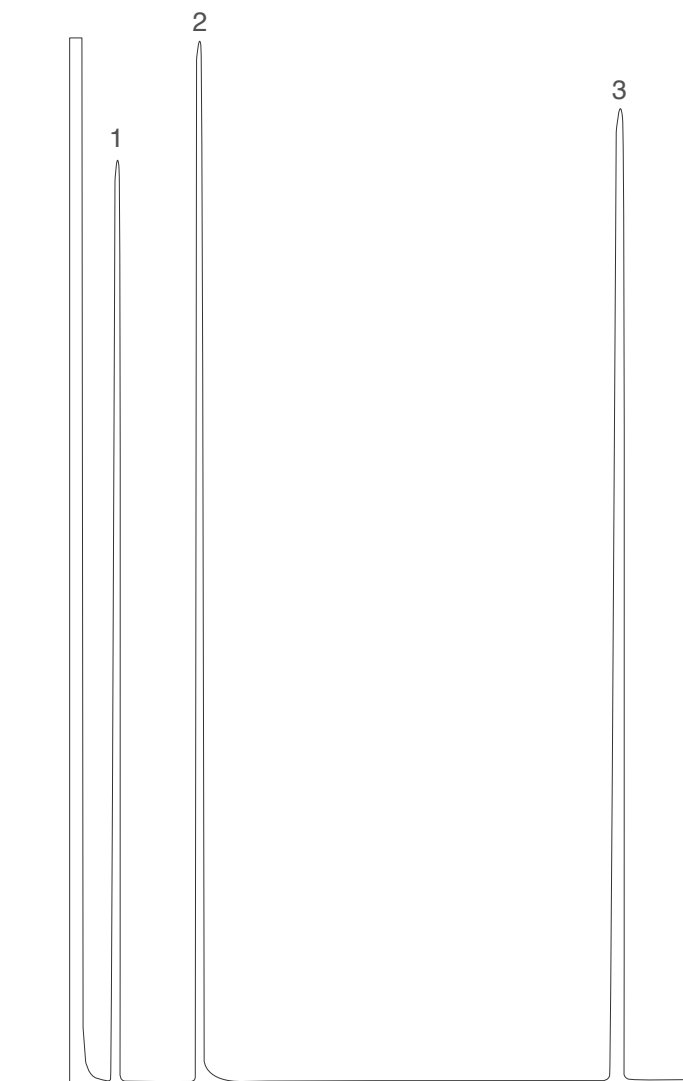

Column: **MEGA-5 - 0.25mm, 0.25 μ m, 25m**

Catalog Code: S-5-025-025-25

Conditions:

Injection: Split 250°C, 1.0 μ L.

Oven Temp: 150°C, 5°C/min, 260°C.

Detector: FID, 280°C

Carrier Gas: Hydrogen, 60 kPa.

Peaks:

1. Amphetamine
2. Ephedrine
3. Caffeine