


# Bergamot Oil



Column: **MEGA-SE54 FAST - 0.10mm, 0.10 $\mu$ m, 5m**

Catalog Code: F-SE54-010-010-5

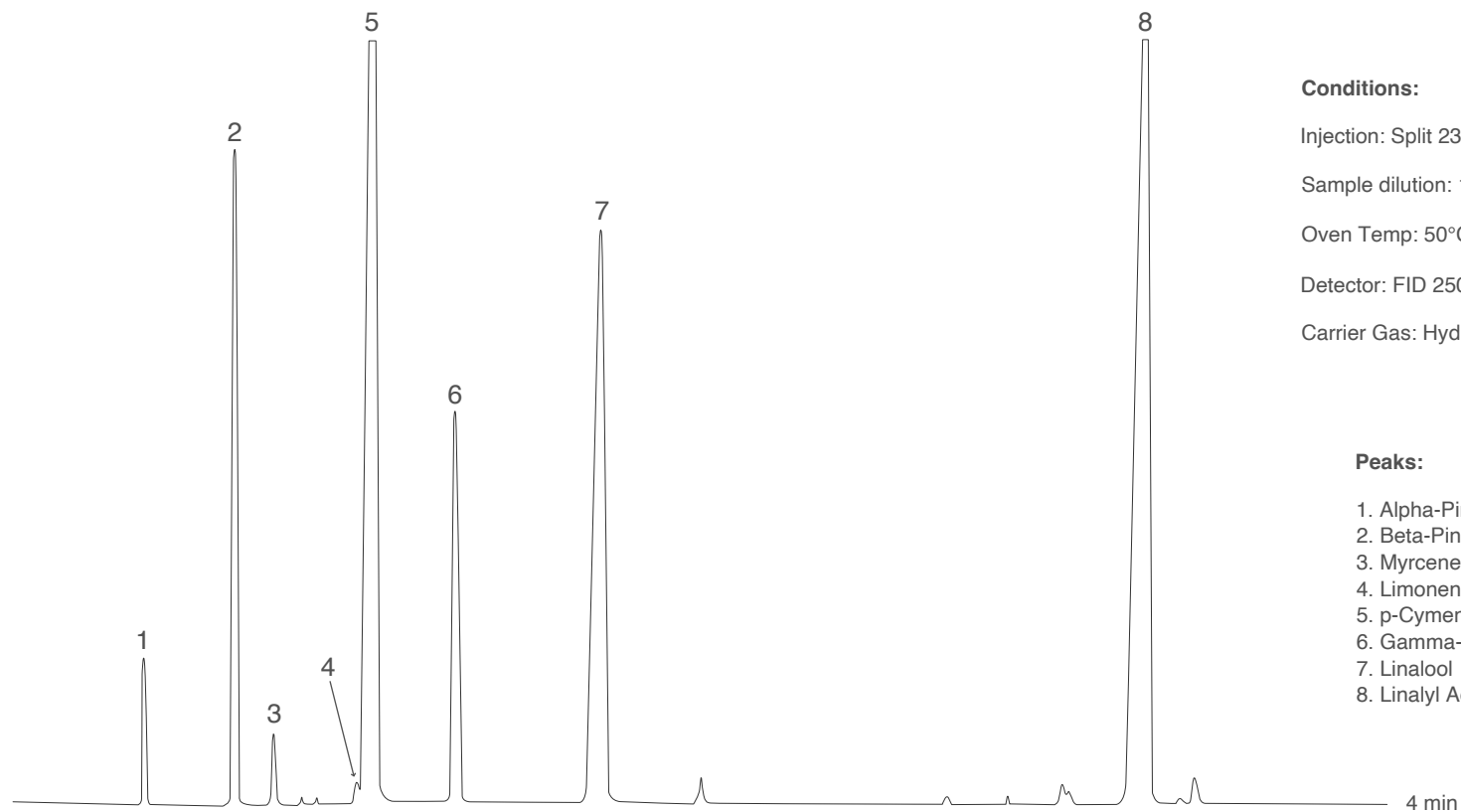

### Conditions:

Injection: Split 230°C, 0.5  $\mu$ L, 1:250 Split Ratio.

Sample dilution: 1% in Cyclohexane.

Oven Temp: 50°C (0.1min), 30°C/min, 250°C.

Detector: FID 250°C

Carrier Gas: Hydrogen, 0.5 mL / min.

### Peaks:

1. Alpha-Pinene
2. Beta-Pinene
3. Myrcene
4. Limonene
5. p-Cymene
6. Gamma-Terpinene
7. Linalool
8. Linalyl Acetate

Acknowledgments: Prof C. Bichi, C. Brunelli et al., Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 Torino, Italy.