


Bergamot Oil

FAST

MEGA[®]
gas chromatography solutions

Column: **MEGA-WAX FAST - 0.10mm, 0.10 μ m, 5m**

Catalog Code: F-WAX-010-010-5

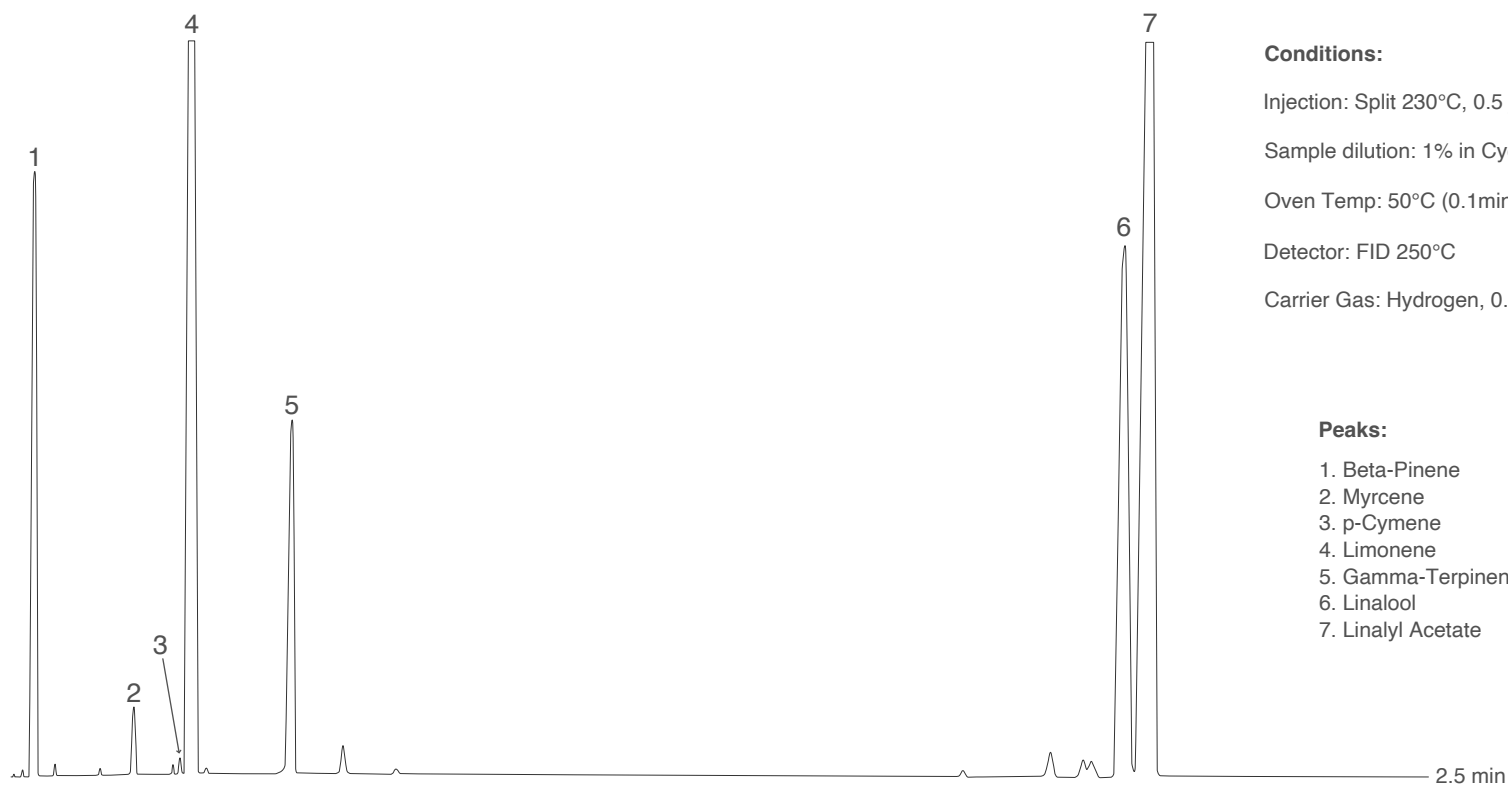

Conditions:

Injection: Split 230°C, 0.5 μ L, 1:250 Split Ratio.

Sample dilution: 1% in Cyclohexane.

Oven Temp: 50°C (0.1min), 30°C/min, 250°C.

Detector: FID 250°C

Carrier Gas: Hydrogen, 0.5 mL / min.

Peaks:

1. Beta-Pinene
2. Myrcene
3. p-Cymene
4. Limonene
5. Gamma-Terpinene
6. Linalool
7. Linalyl Acetate

Acknowledgments: Prof C. Bicchi, C. Brunelli et al., Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 Torino, Italy.