


# Butter Triglycerides (T24 - T56)

STANDARD

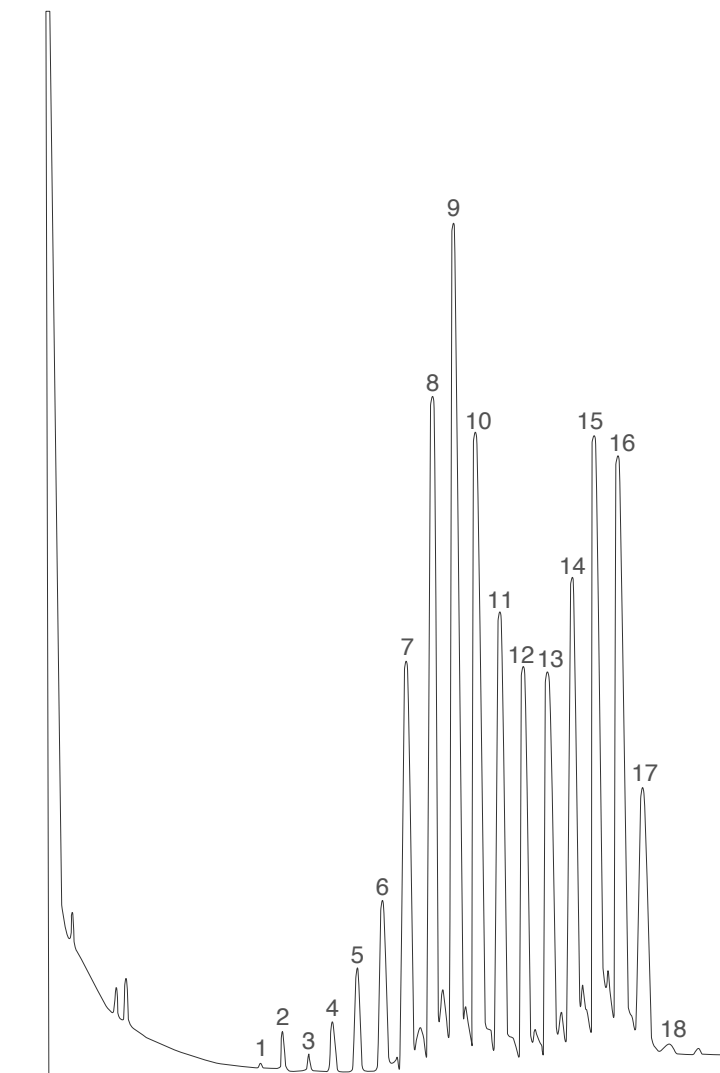

Column: **MEGA-1 - 0.32mm, 0.10 $\mu$ m, 3m + 1m built-in Retention Gap**

Catalog Code: GAP-1-032-010-3-1

**Conditions:**

Injection: On-Column.

Oven Temp: 60°C, 15°C /min, 340°C.

Detector: FID, 350°C.

Carrier Gas: Hydrogen 30 kPa.

**Peaks:**

- |                |         |
|----------------|---------|
| 1. T24         | 11. T42 |
| 2. Cholesterol | 12. T44 |
| 3. T26         | 13. T46 |
| 4. T28         | 14. T48 |
| 5. T30         | 15. T50 |
| 6. T32         | 16. T52 |
| 7. T34         | 17. T54 |
| 8. T36         | 18. T56 |
| 9. T38         |         |
| 10. T40        |         |

T = Total number of carbons