


Column: **MEGA-WAX FAST - 0.10mm, 0.10 μ m, 5m**
Catalog Code: F-WAX-010-010-5

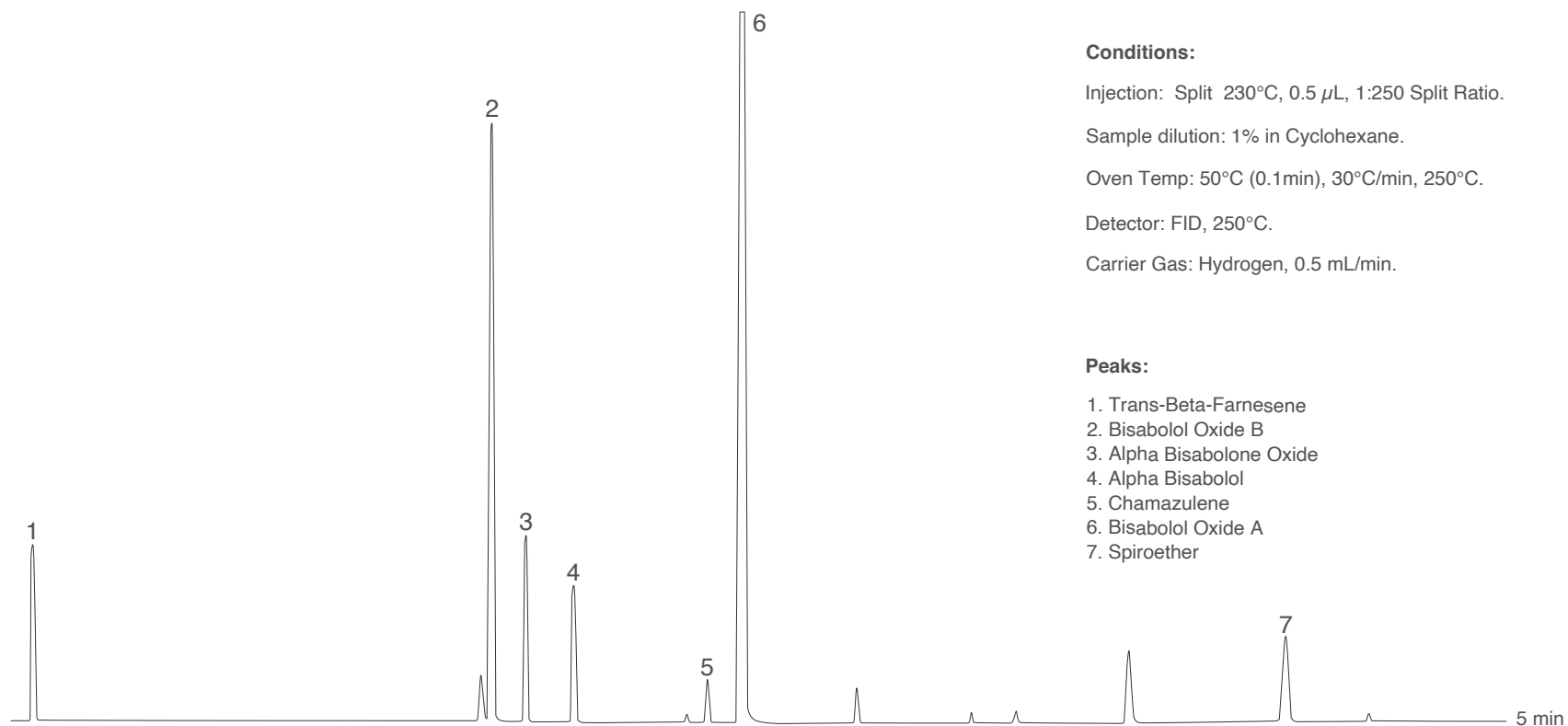

Conditions:

Injection: Split 230°C, 0.5 μ L, 1:250 Split Ratio.

Sample dilution: 1% in Cyclohexane.

Oven Temp: 50°C (0.1min), 30°C/min, 250°C.

Detector: FID, 250°C.

Carrier Gas: Hydrogen, 0.5 mL/min.

Peaks:

1. Trans-Beta-Farnesene
2. Bisabolol Oxide B
3. Alpha Bisabolone Oxide
4. Alpha Bisabolol
5. Chamazulene
6. Bisabolol Oxide A
7. Spiroether

Acknowledgement: Prof C. Bicchi, C. Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.