

Alpha-Pinene (+/-), Linalool (+/-) and Menthol (+/-)


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DEX DET Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DET-B-025-025-25

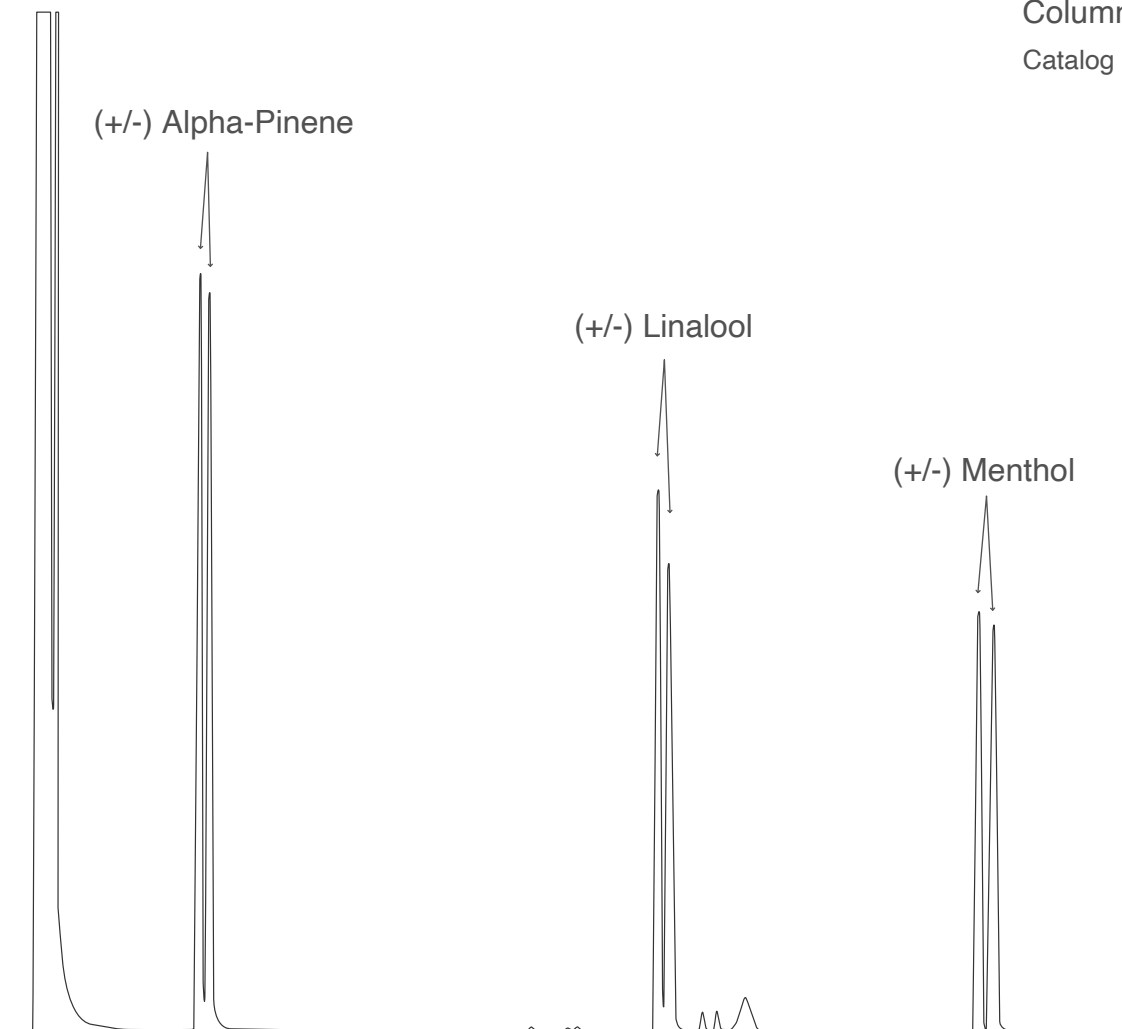

Conditions:

Injection: Split 250°C, 0.5 μ L, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.

Oven Temp.: 75°C, 1°C/min, 160°C.

Detector: FID, 250°C.

Carrier Gas: Hydrogen, 75 kPa.

CHIRAL DEX