

Delta-Lactones

Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DAC-B-025-025-25

Peaks:


1. (+/-) Delta-Octalactone
2. (+/-) Delta-Nonalactone
3. (+/-) Delta-Decalactone
4. (+/-) Delta-Undecalactone
5. (+/-) Delta-Dodecalactone

Conditions:

Injection: Split 250°C, 0,5 μ L, 1:50 Split Ratio.

Sample dilution: 1mg/mL each in n-Hexane.

Contact us for additional information.

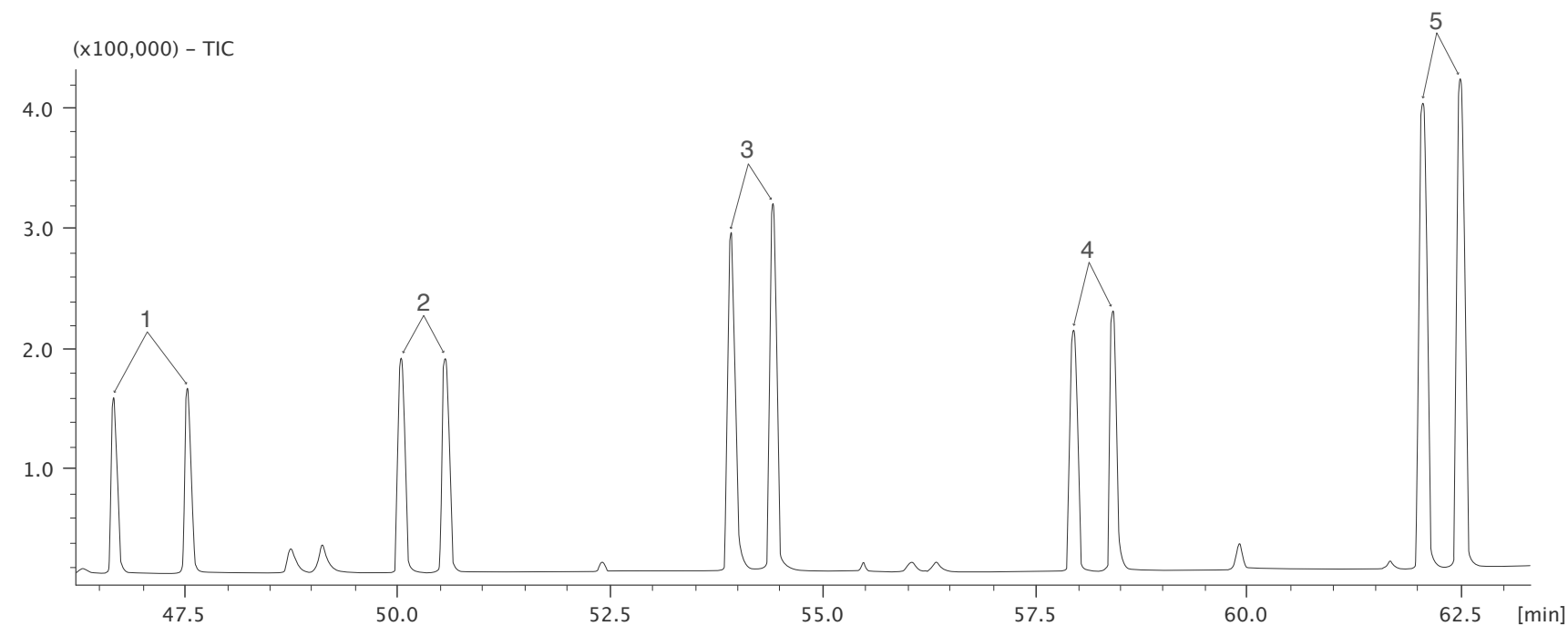

Acknowledgement: Prof C. Bicchi, C. Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.