


Tryglicerides - Hazelnut Oil

CUSTOM
DEDICATED

MEGA®
gas chromatography solutions

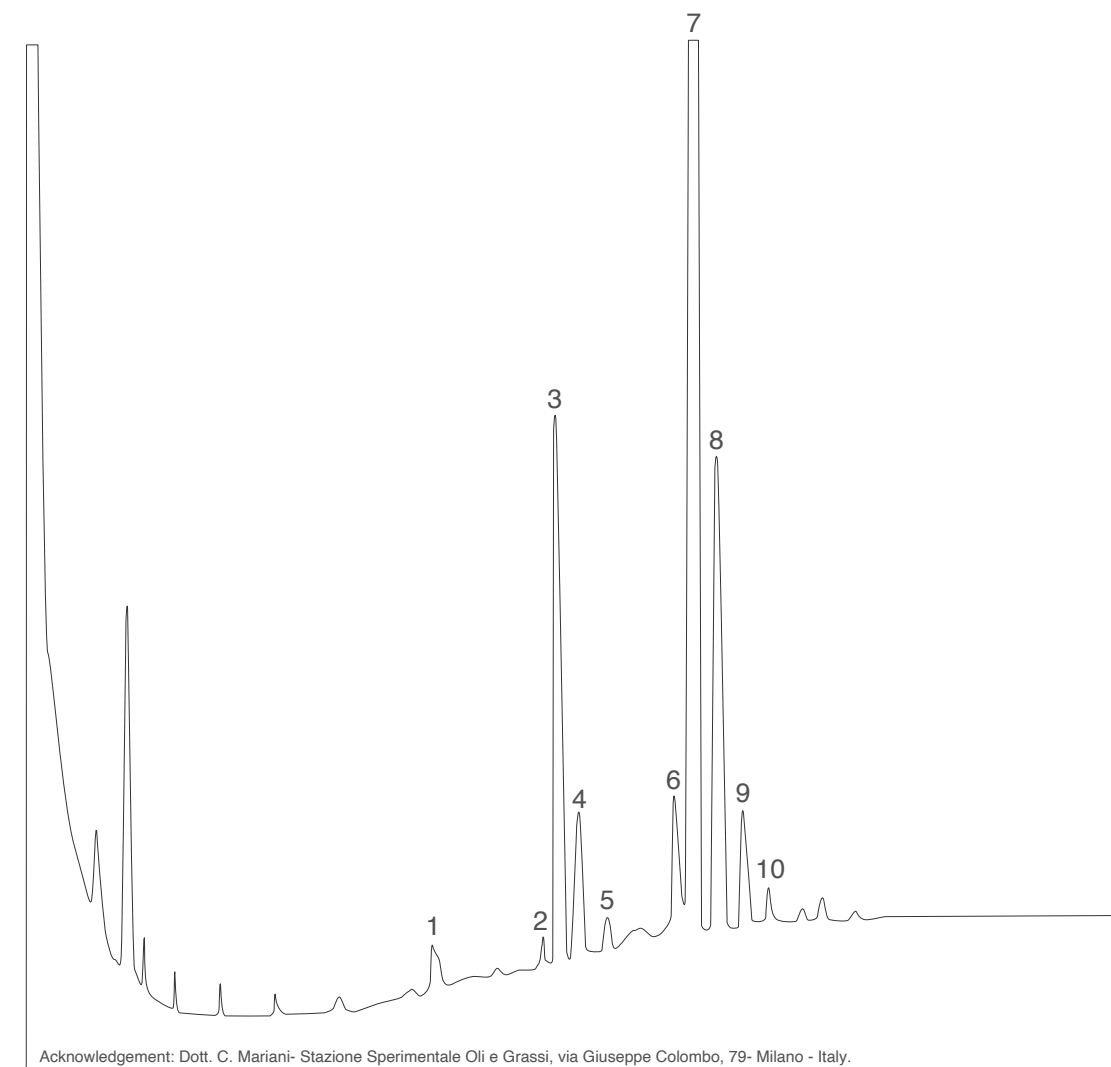

Column: **MEGA-LAP - 0.32mm, 0.10 μ m, 25m**

Catalog Code: C-LAP-032-010-25

Conditions:

Injection: Split 300°C, 1.0 μ L, 1:80 Split Ratio.

Detector: FID, 300°C.

Oven Temp.: 120°C, 5°C/min, 330°C.

Carrier Gas: Hydrogen, 2 mL/min.

Peaks:

1. PPO
2. POS
3. POO
4. PLO+OOPo
5. PLL+PoOI
6. SOO
7. OOO
8. OOL
9. OLL
10. LLL