

Tryglicerides - Olive Oil

CUSTOM
DEDICATED

MEGA®
gas chromatography solutions

Column: **MEGA-LAP - 0.32mm, 0.10µm, 25m**

Catalog Code: C-LAP-032-010-25

Conditions:

Injection: Split 300°C, 1.0 µL, 1:80 Split Ratio.


Detector: FID, 300°C.

Oven Temp.: 120°C, 5°C/min, 330°C.

Carrier Gas: Hydrogen, 2 mL/min.

Peaks:

1. PPO
2. POS
3. POO
4. PLO+OOPo
5. PLL+PoOI
6. SSO
7. SSO
8. OOO
9. OOL
10. OLL
11. LLL
12. AOO
13. GoOO

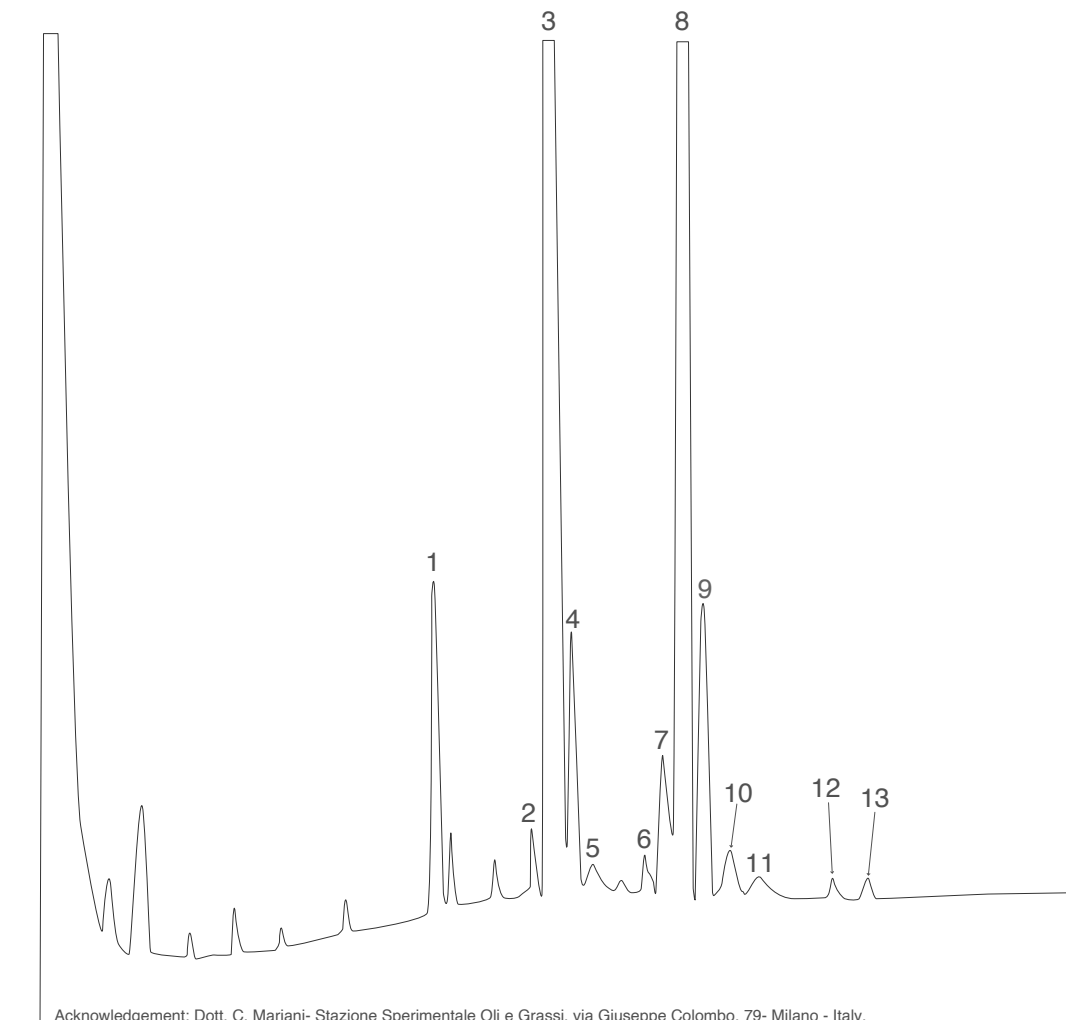

Acknowledgement: Dott. C. Mariani- Stazione Sperimentale Oli e Grassi, via Giuseppe Colombo, 79- Milano - Italy.