

Vodka (zoom first part of the chromatogram)

View the related entire Vodka chromatogram

Column: **MEGA-WAX Spirit - 0.32mm, 0.50 μ m, 50m**

Catalog Code: C-WAXSP-032-050-50


Conditions:

Injection: Split 250°C, 1.0 μ L, 50 mL/min Split Ratio.

Detector: FID, 250°C.

Oven Temp.: 60°C (2min), 10°C/min, 250°C.

Carrier Gas: Nitrogen, 0.8 mL/min.

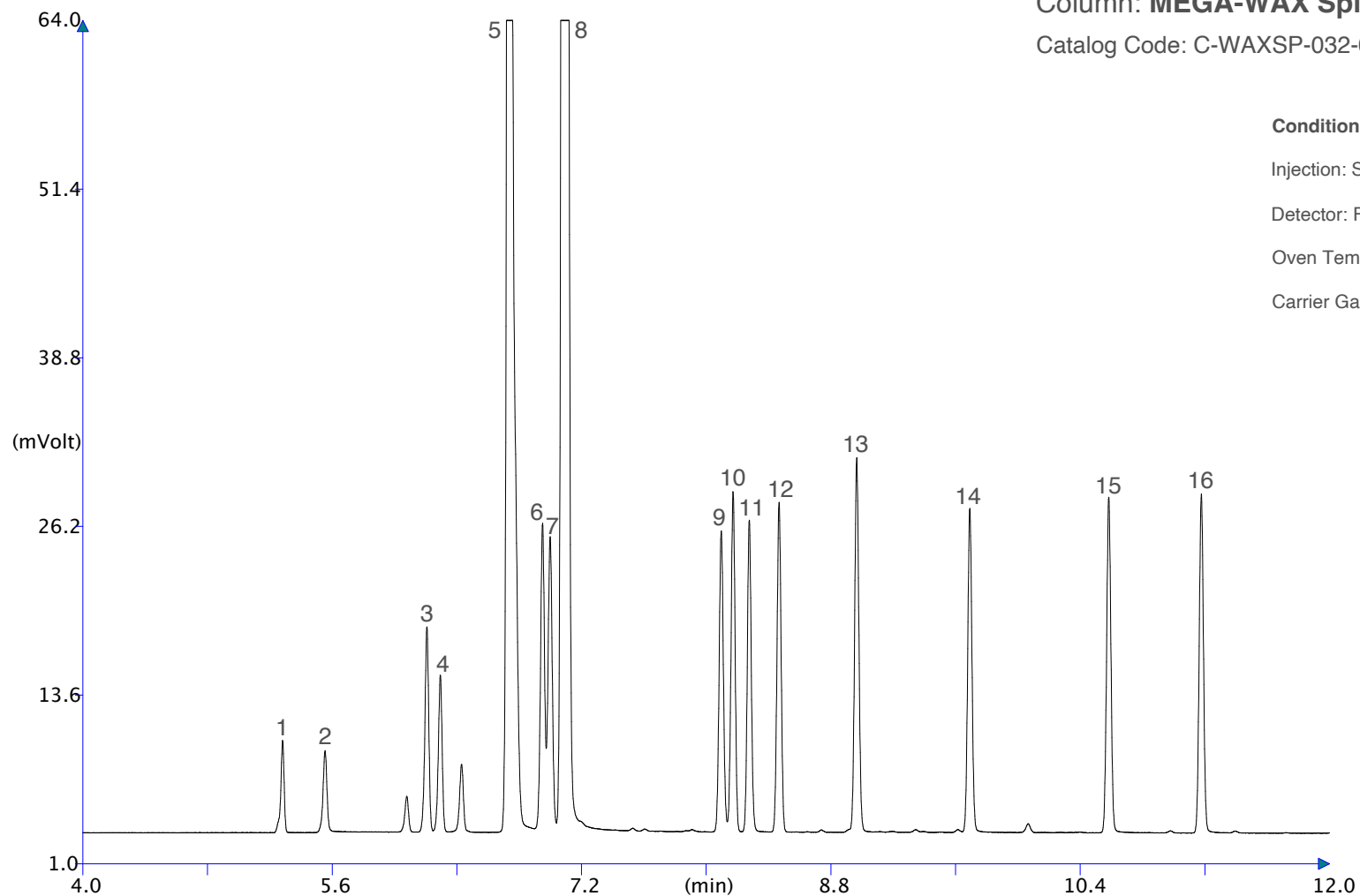

Peaks:

1. Diethyl ether
2. Acetaldehyde
3. Acetone
4. Methylacetate
5. Methanol
6. MEK
7. 2-propanol
8. Ethanol
9. Isobutylacetate
10. 2-butanol
11. 1-propanol
12. Ethylbutyrate
13. Crotonic aldehyde
14. Isobutanol
15. 1-butanol
16. Isopentanol