


Volatile Acidic Compounds

CUSTOM
DEDICATED

MEGA[®]
gas chromatography solutions

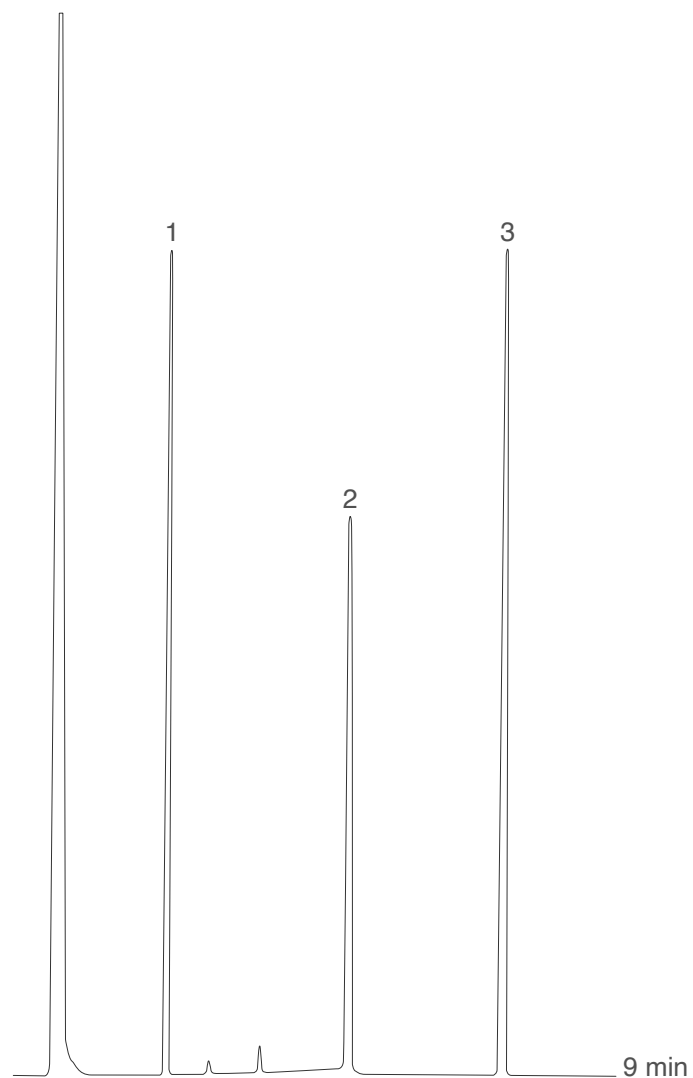

Column: **MEGA-ACID (FFAP) - 0.53mm, 1.00 μ m, 10m**

Catalog Code: C-ACID-053-100-10

Conditions:

Injection: Split 250°C, 1.0 μ L.

Detector: FID, 250°C.

Oven Temp.: 80°C, 2°C/min, 120°C.

Carrier Gas: Helium, 20 kPa.

Peaks:

1. Acetic anhydride
2. Acetic acid
3. Isobutyric acid