

# 1,2-Epoxydodecane


CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Chiral separation on MEGA-DEX columns

Column: **MEGA-DEX ASX 2-92 - 0.25mm, 0.25 $\mu$ m, 30m**

Catalog Code: DEX-ASX2-92-025-025-30

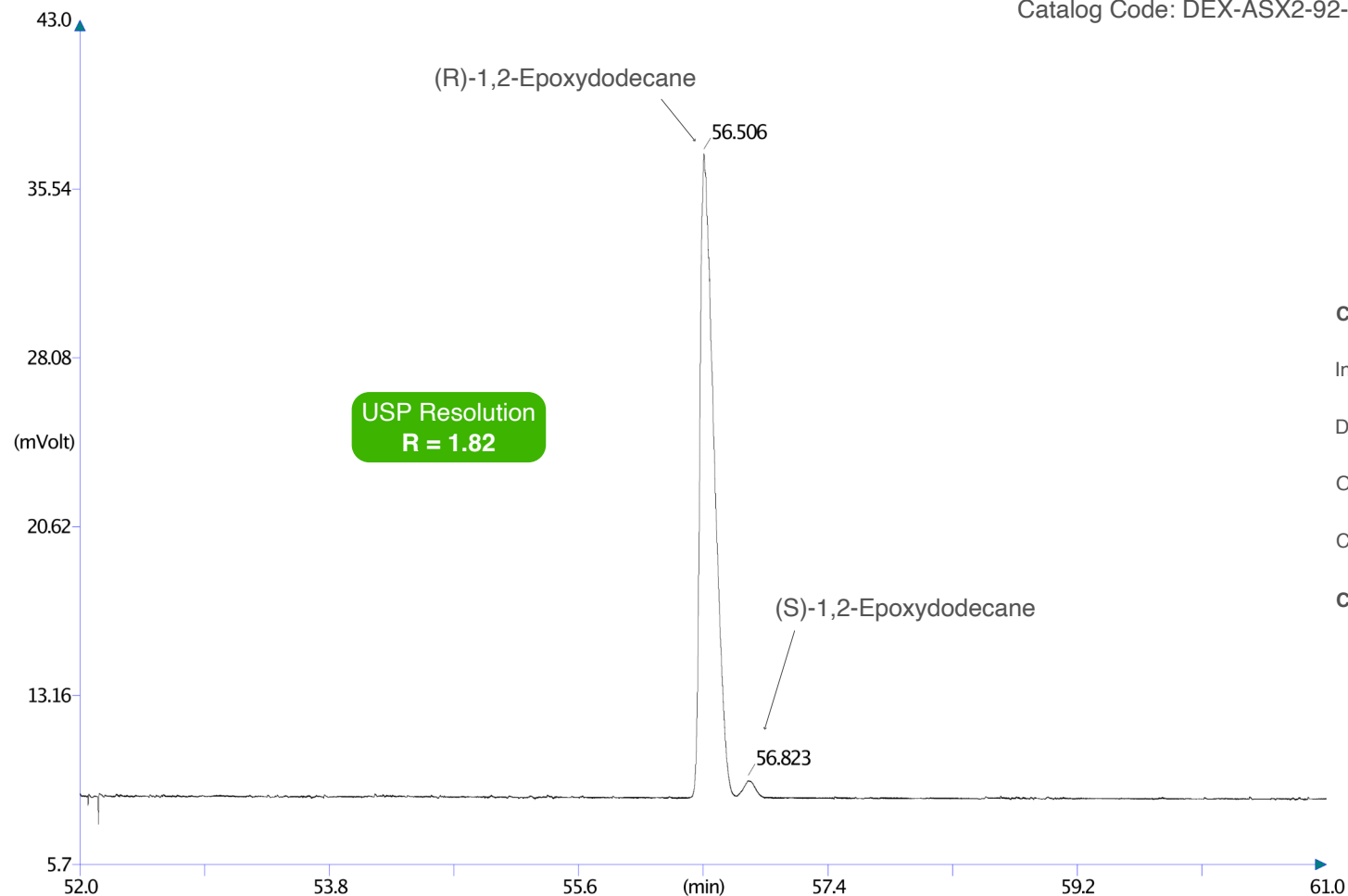

#### Conditions:

Injection: Split, 250°C.

Detector: FID, 250°C.

Oven Temp.: 50°C, 1.5°C/min.

Carrier Gas: Hydrogen, 90 kPa.

Contact us for additional information.