

(R/S) 1,2-Propanediol


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 30m**

Catalog Code: DEX-DMT-B-025-025-30

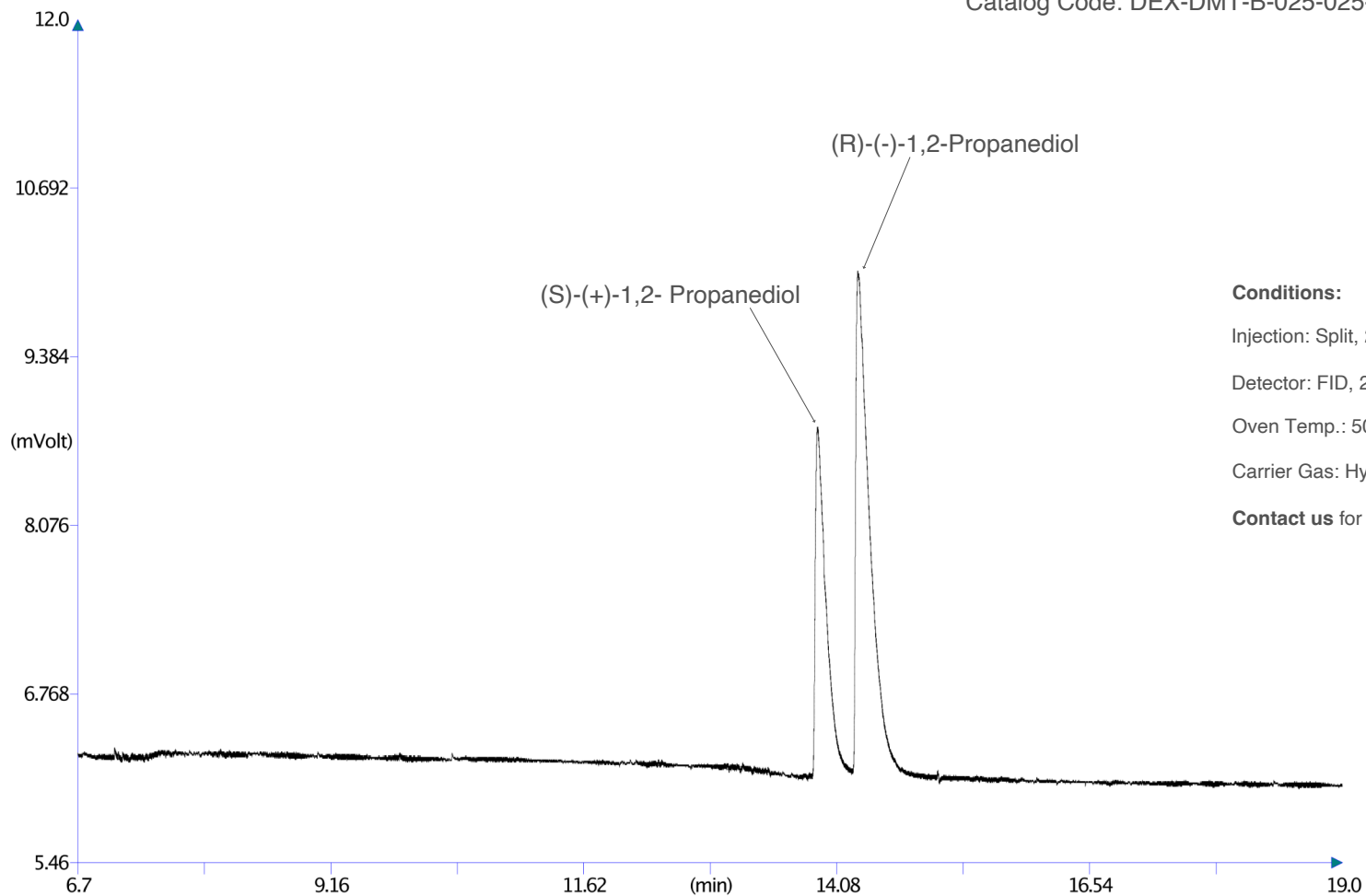

Conditions:

Injection: Split, 250°C.

Detector: FID, 250°C.

Oven Temp.: 50°C, 1.5°C/min.

Carrier Gas: Hydrogen, 80 kPa.

Contact us for additional information.