

2-oxiranecarboxylic acid, 3-propyl, methyl ester (+/-)


CHIRAL DEX

MEGA[®]
gas chromatography solutions

Chiral separation on MEGA-DEX columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

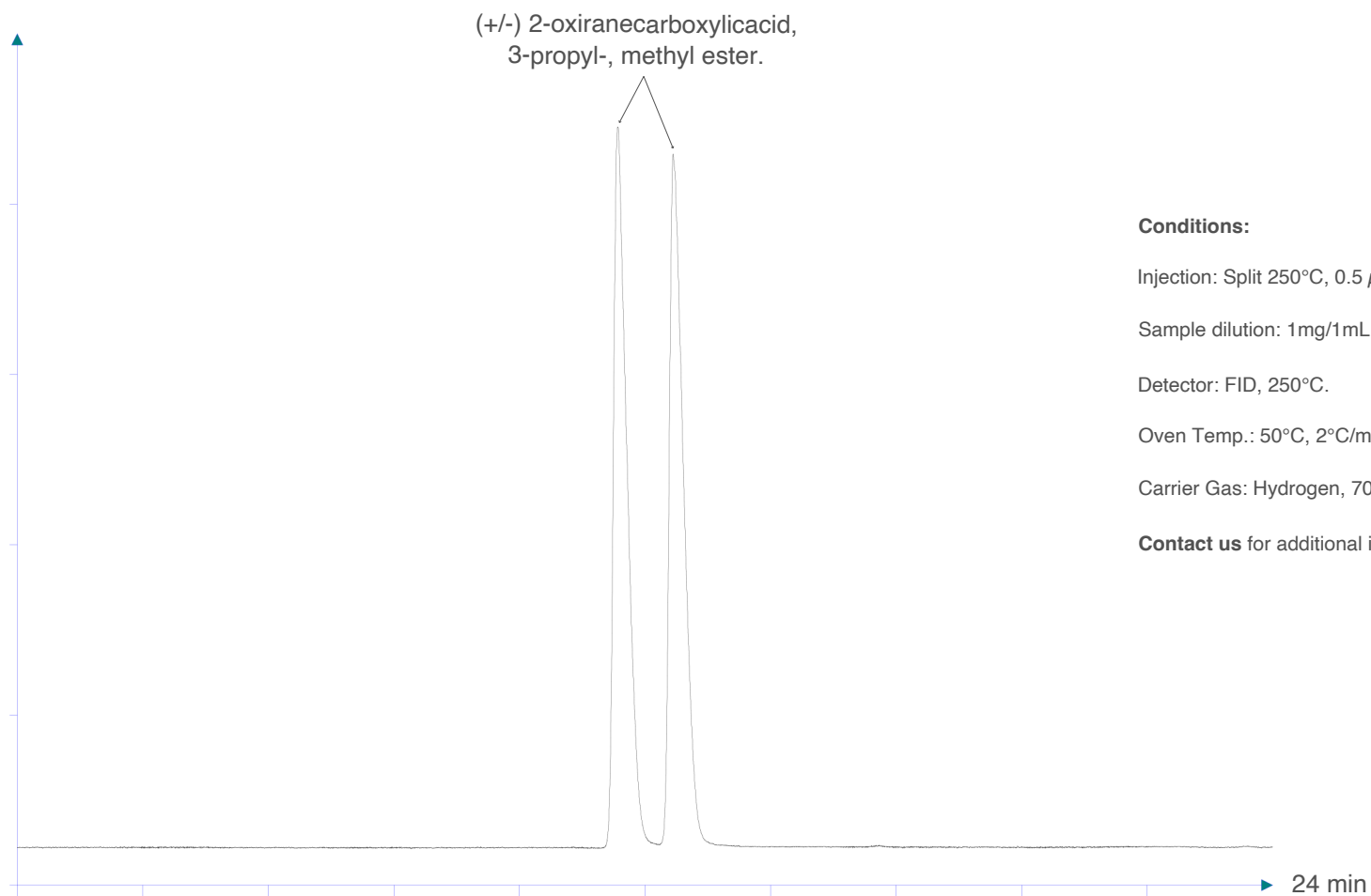

Conditions:

Injection: Split 250°C, 0.5 μ L, 90 mL/min Split Flow.

Sample dilution: 1 mg/1 mL in Methylene Chloride.

Detector: FID, 250°C.

Oven Temp.: 50°C, 2°C/min, 200°C.

Carrier Gas: Hydrogen, 70 kPa.

Contact us for additional information.