

# Alpha-Copaene (+/-)

Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Column: **MEGA-DEX DMP Beta - 0.25mm, 0.25 $\mu$ m, 25m**

Catalog Code: DEX-DMP-B-025-025-25

**Conditions:**

Injection: Split, 250°C.


Detector: MS.

Oven Temp.: 50°C, 2°C/min.

Carrier Gas: Helium, 1mL/min.

**Contact us** for additional information.

All our MEGA-DEX chiral columns work perfectly with GC-MS systems.

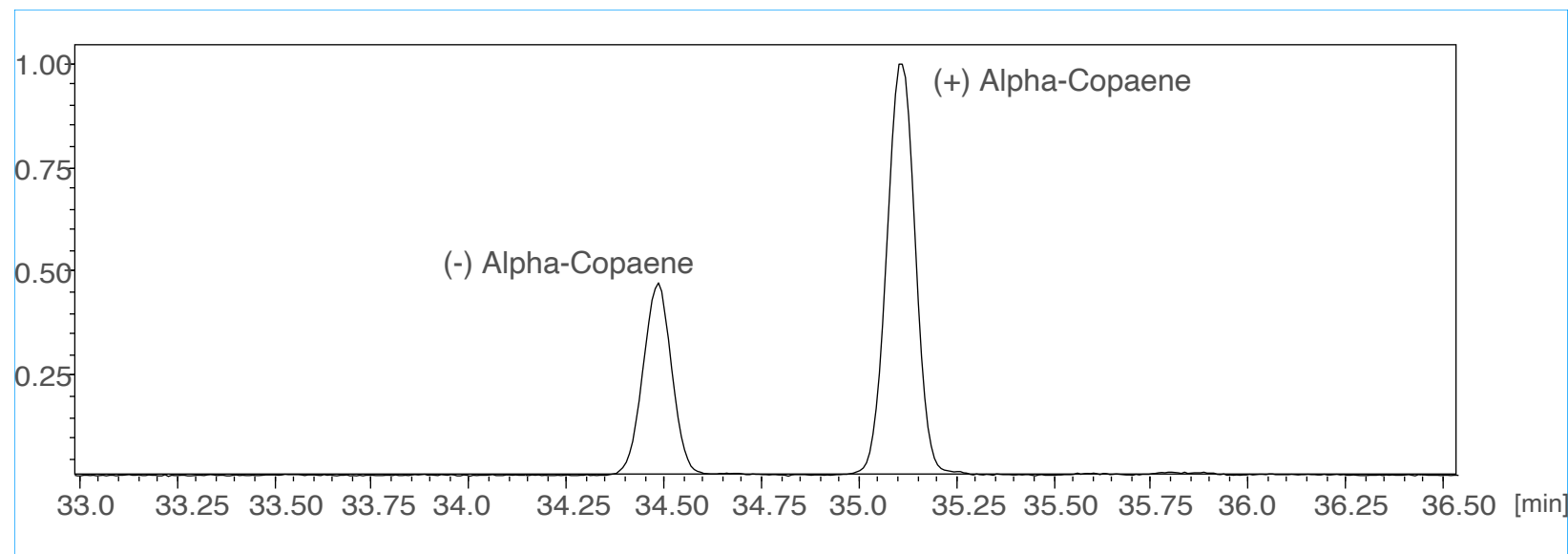

Acknowledgement: many thank to Barbara Sgorbini, Cecilia Cagliero, Chiara Cordero and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, via P. Giuria,9, Torino, Italy.