


Benzene - 1,2-Dichloroethane - MEGA-VOC 1 & 2 comparison



Column: **MEGA-VOC 2 - 0.32mm, 5.00 μ m, 60m**

Catalog Code: C-VOC2-032-500-60

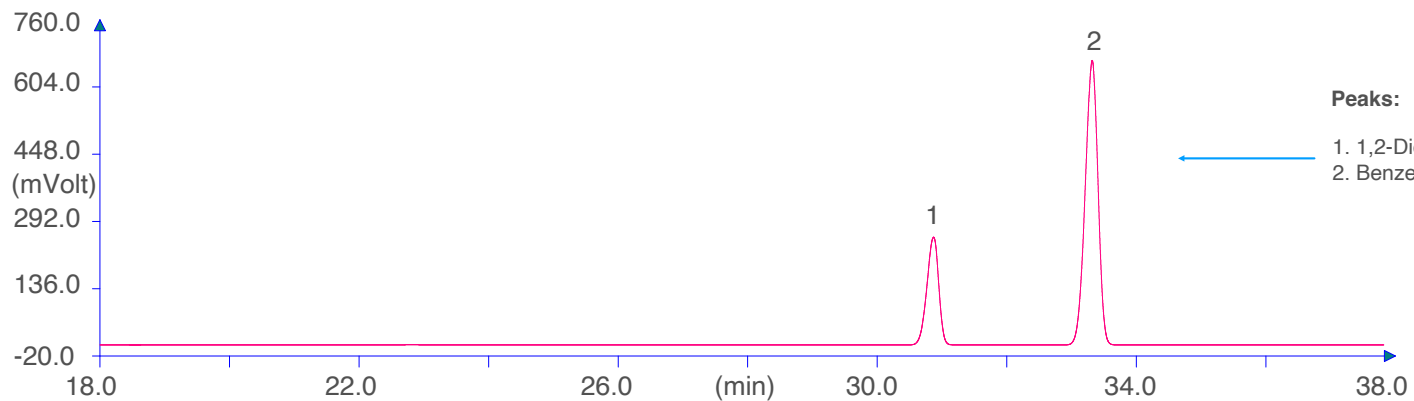

Peaks:

- 1. 1,2-Dichloroethane (31.05 min)
- 2. Benzene (33.50 min)

Conditions:

Injection: Split, 250°C.


Detector: FID, 250°C.

Oven Program: 35°C (5min), 1°C/min.

Carrier Gas: Helium, 120kPa.

Column: **MEGA-VOC 1 - 0.32mm, 3.00 μ m, 60m**

Catalog Code: C-VOC1-032-300-60

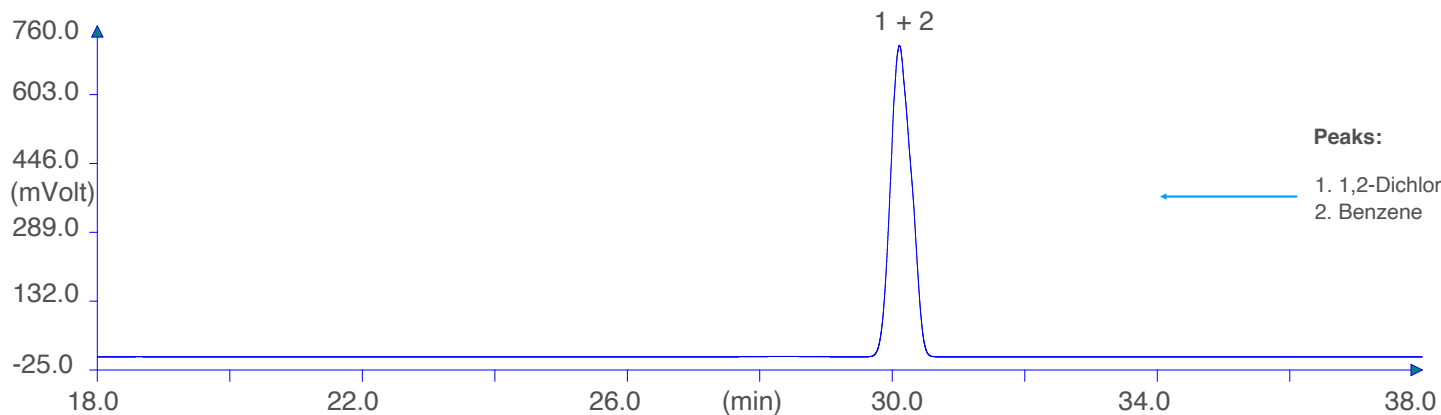

Peaks:

- 1. 1,2-Dichloroethane
- 2. Benzene