

## Beta-Citronellol (+/-)


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25 $\mu$ m, 25m**

Catalog Code: DEX-DAC-B-025-025-25

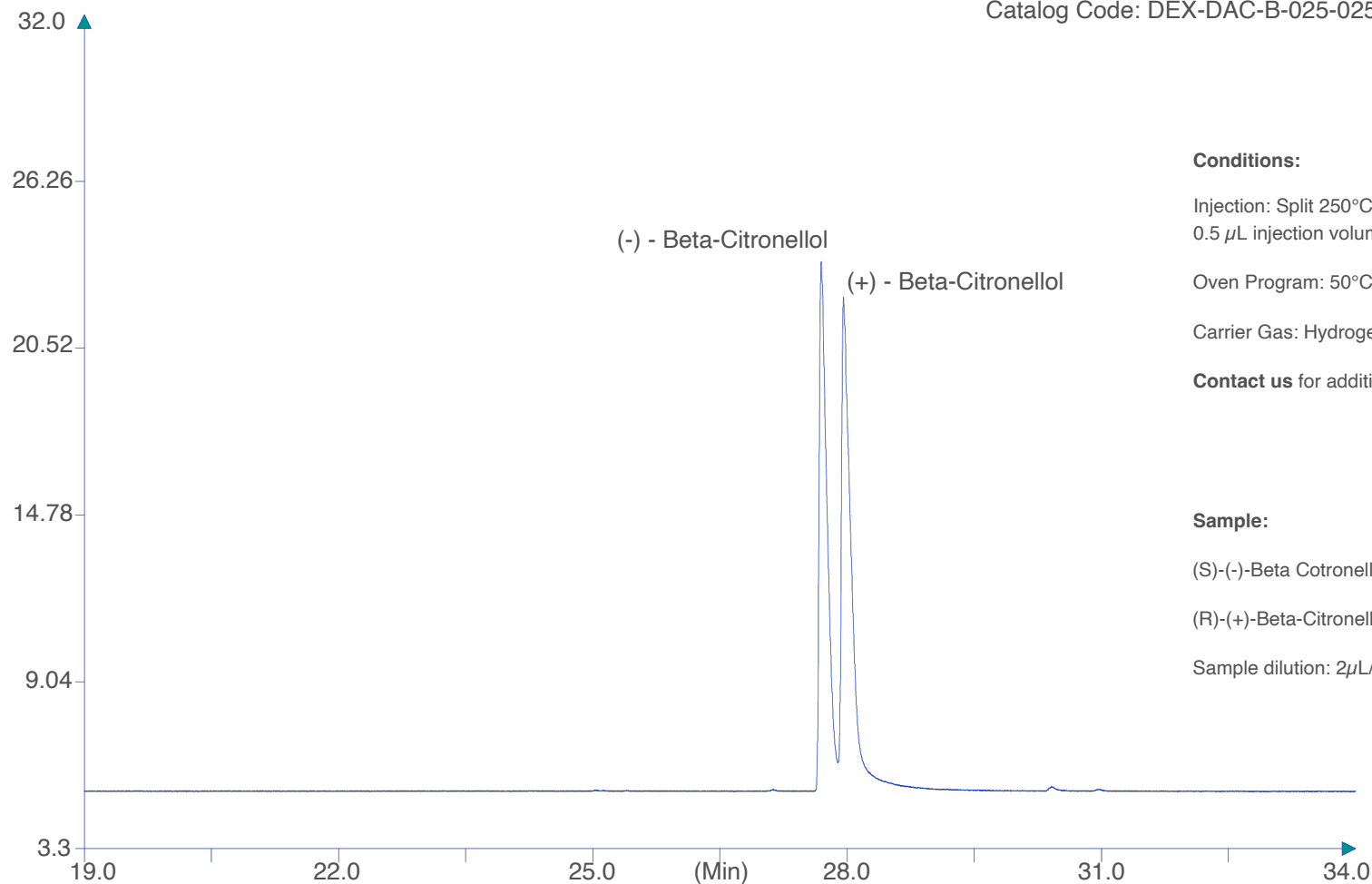

### Conditions:

Injection: Split 250°C, 100mL/min Split flow,  
0.5  $\mu$ L injection volume.

Oven Program: 50°C, 2°C/min.

Carrier Gas: Hydrogen, 80 kPa.

**Contact us** for additional information.

### Sample:

(S)-(-)-Beta Citronellol (CAS# 7540-51-4).

(R)-(+)-Beta-Citronellol (CAS# 1117-61-9).

Sample dilution: 2 $\mu$ L/1mL in n-Hexane each.