


Food Packaging Residual Solvents



On MEGA-VOC 1 & MEGA-VOC 2 dual system

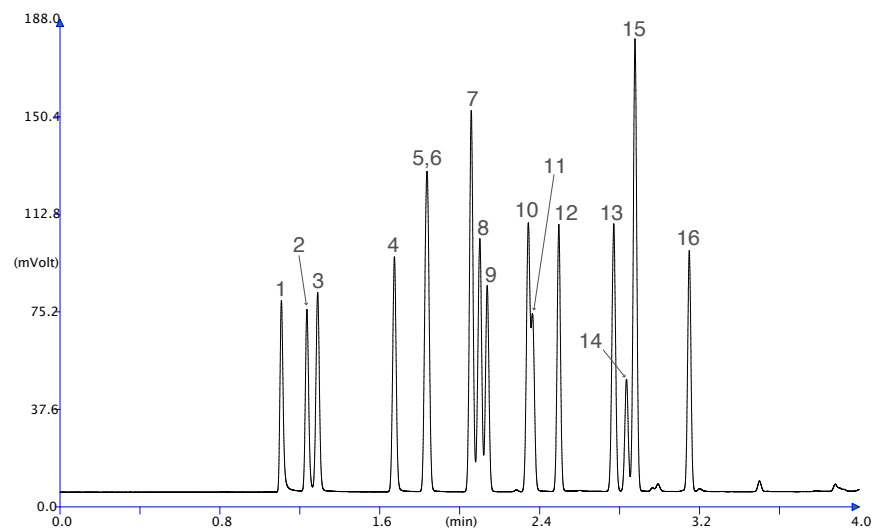

Column: **MEGA-VOC 1 FAST - 0.20mm, 1.10 μ m, 10m**

Catalog Code: F-VOC1-020-110-10


Conditions:

Injection: Split 250°C, 0,2 μ L (with micro-syringe, 150mL/min Split Flow.

Oven Temp.: 35°C (0.4 min), 35°C/ min, 200°C.

Detector: FID, 250°C.

Carrier Gas: Helium, 60 kPa.

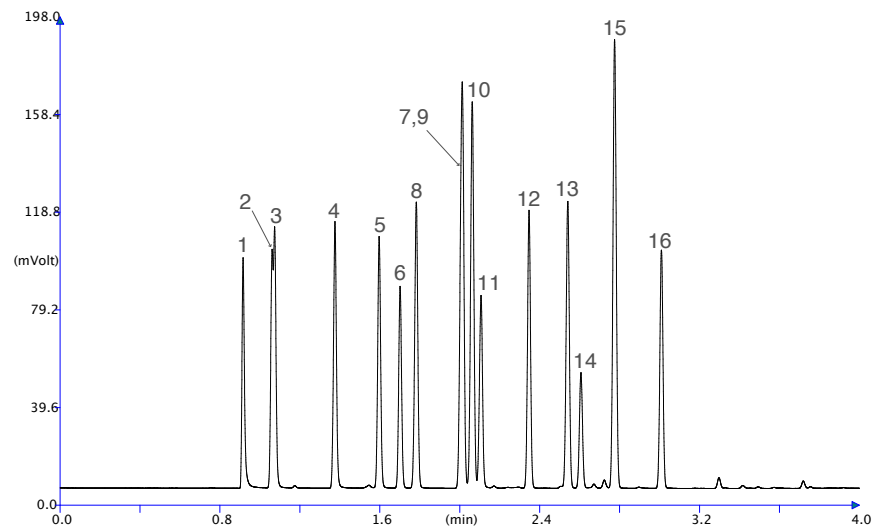

Column: **MEGA-VOC 2 FAST - 0.20mm, 1.00 μ m, 10m**

Catalog Code: F-VOC2-020-100-10

Peaks:

- | | |
|------------------------|----------------------------|
| 1. Ethanol | 9. Isopropyl acetate |
| 2. Acetone | 10. 1-butanol |
| 3. 2-propanol | 11. Methoxy propanol |
| 4. 1-propanol | 12. Propyl acetate |
| 5. Methyl-Ethyl-Ketone | 13. Isobutyl-Methyl-Ketone |
| 6. Ethyl acetate | 14. Ethoxy propanol |
| 7. Cyclohexane | 15. Toluene |
| 8. Isobutanol | 16. Butyl acetate |