

Hexobarbital (Pharmaceuticals)


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DEX DMP Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMP-B-025-025-25

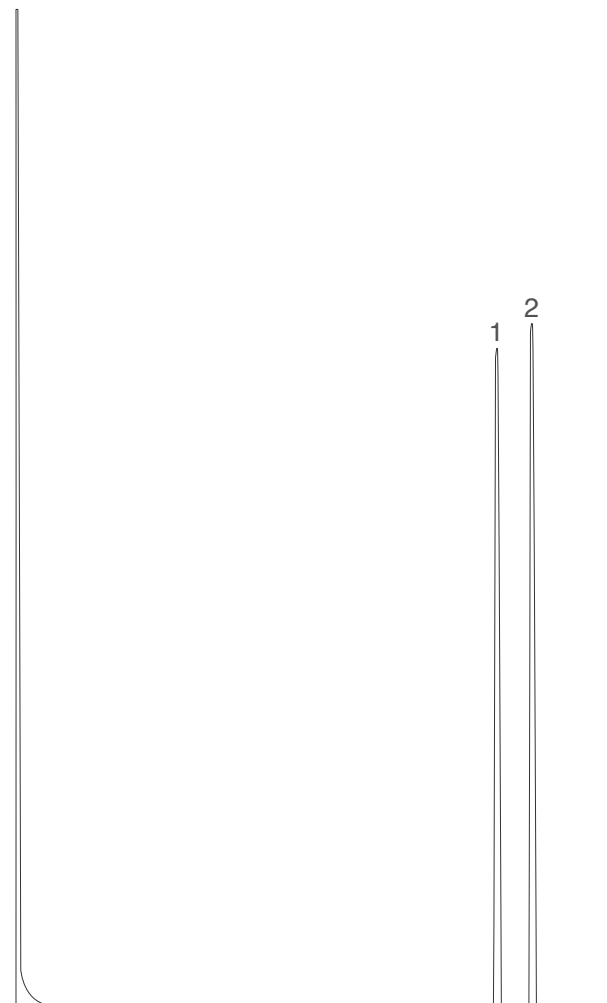

Conditions:

Injection: Split 250°C, 0.5 μ L, 1:50 Split Ratio.

Sample dilution: 1 mg/mL in n-Hexane.

Detector: FID, 250°C.

Oven Temp.: 100°C, 2.5°C/min, 170°C.

Carrier Gas: Hydrogen, 70 kPa.

Peaks:

1. (+) Hexobarbital
2. (-) Hexobarbital