

# Linalool, Linalool Oxyde

Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Column: MEGA-DEX DMT Beta - 0.25mm, 0.25 $\mu$ m, 25m


Catalog Code: DEX-DMT-B-025-025-25

## Conditions:

Injection: Split 250°C, 0.5  $\mu$ L, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.

Contact us for additional information.

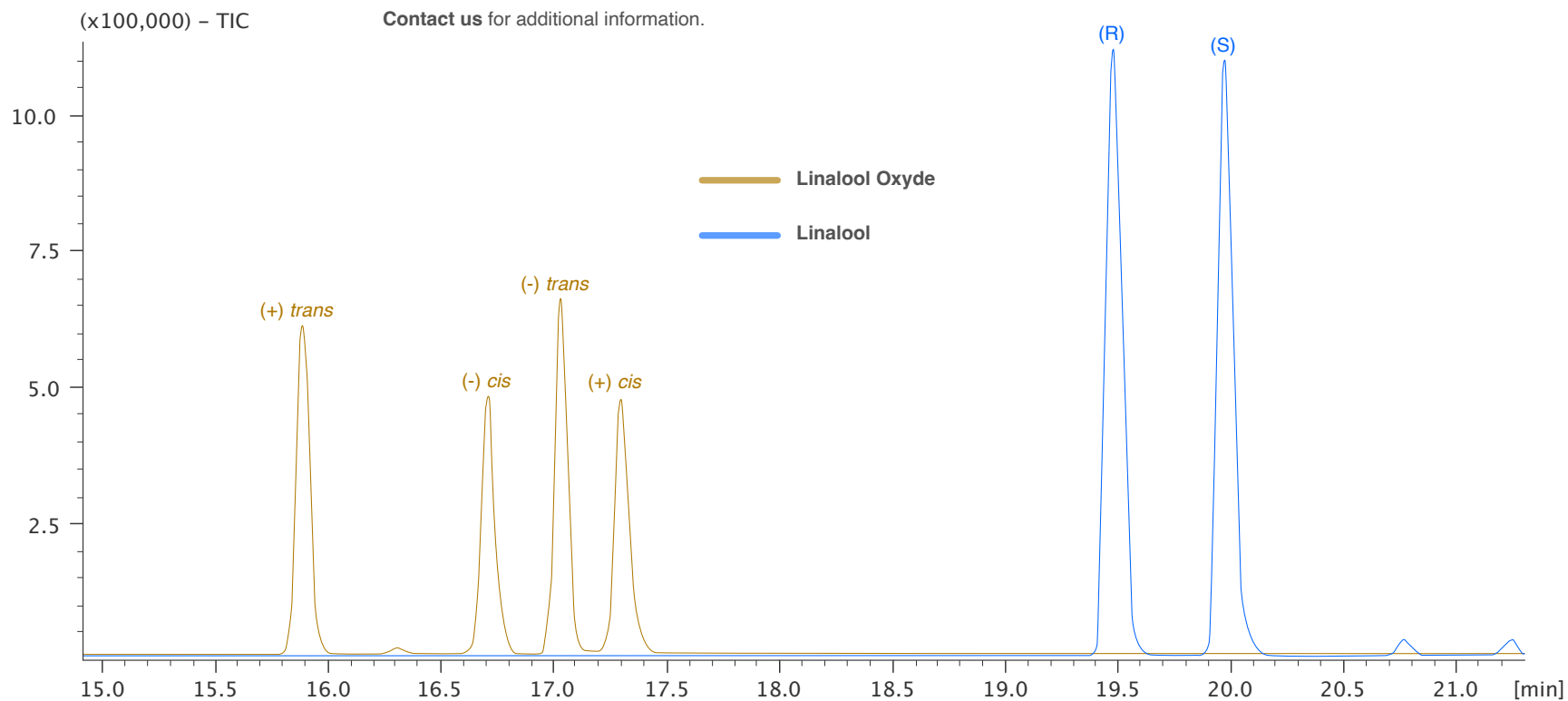

Acknowledgement: Prof C. Bicchi, C. Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.