

# (Iso)Menthol, Neomenthol and Neoisomenthol

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions

Chiral separation on MEGA-DEX columns

Column: **MEGA-DEX DAC Beta - 0.25mm, 0.25 $\mu$ m, 25m**


Catalog Code: DEX-DAC-B-025-025-25

### Conditions:

Injection: Split 250°C, 0.5  $\mu$ L, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.

Contact us for additional information.

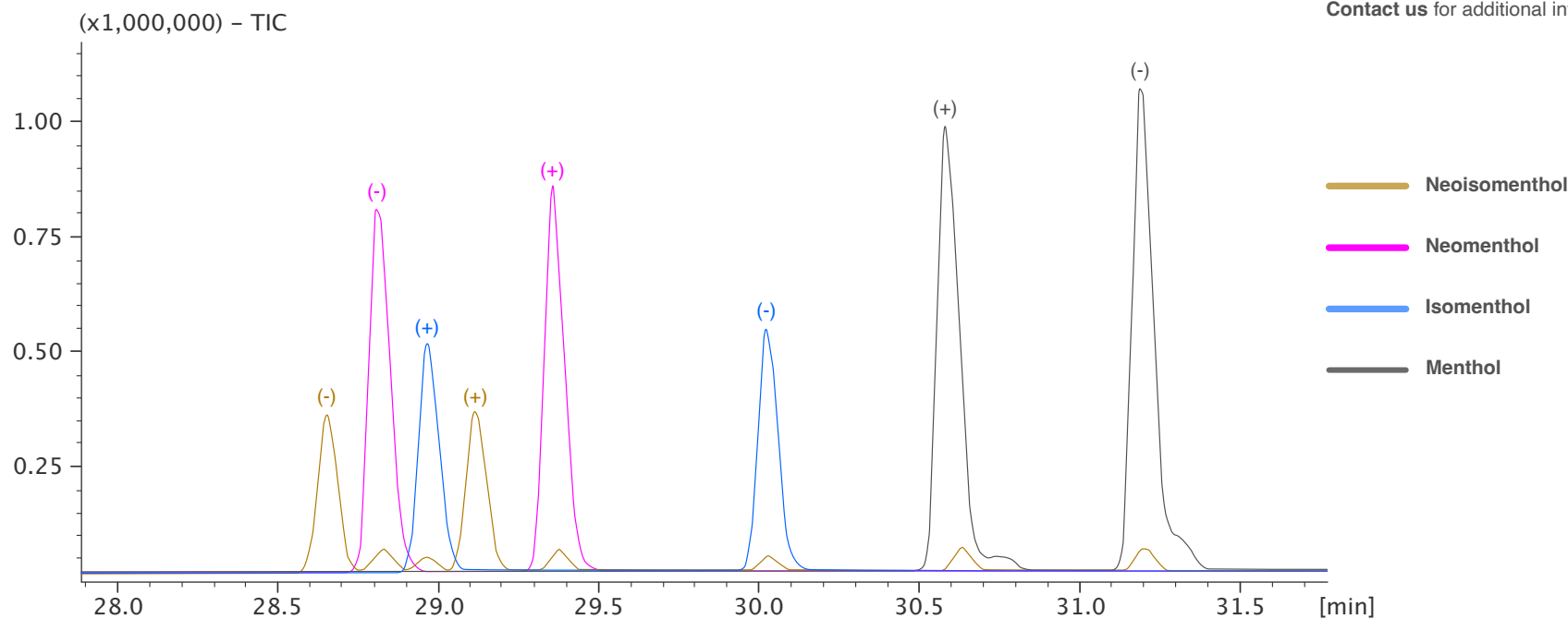

Acknowledgement: Prof C. Bicchi, C. Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.