


Nitrosamines

STANDARD

MEGA[®]
gas chromatography solutions

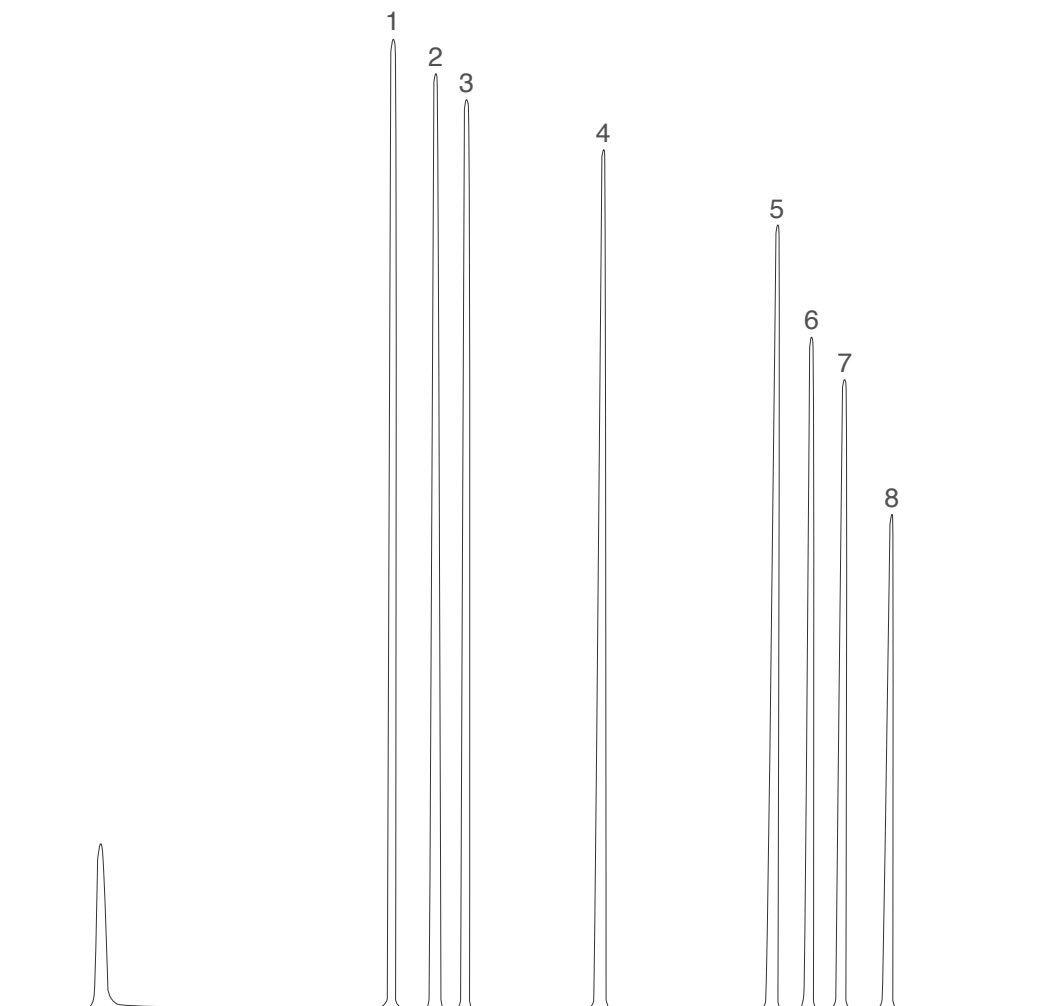

Column: **MEGA-WAX - 0.32mm, 1.20 μ m, 25m**

Catalog Code: S-WAX-032-120-25

Conditions:

Injection: Spitless, 200°C.

Detector: NPD.

Oven Temp.: 60°C (1min), 0.5°C/min, 130°C (2min),
25°C/min, 160°C.

Carrier Gas: Helium, 60 kPa.

Peaks:

1. N-nitrosodimethylamine
2. N-nitrosomethylethylamine
3. N-nitrosodiethylamine
4. N-nitrosodi-n-propylamine
5. N-nitrosodi-n-butylamine
6. N-nitrosopiperidine
7. N-nitrosopirrolidine
8. N-nitrosomorpholine