

Column: **MEGA-1701 FAST - 0.10mm, 0.10 μ m, 5m**

Catalog Code: F-1701-010-010-5

Conditions:

Injection: Split 230°C, 0.5 μ L, 1:200 Split Ratio.

Sample: 1% in Cyclohexane.


Detector: FID 250°C.

Oven Temp.: 50 °C(0.1min), 40°C/min, 250°C.

Carrier Gas: Hydrogen, 0,5mL/min.

Peaks:

- | | |
|---------------------|---------------------------------|
| 1. Alpha-HCH | 10. / |
| 2. Gamma HCH | 11. Fenitrothion |
| 3. Heptachlor | 12. Chlordane- <i>cis+trans</i> |
| 4. Chlorothalonil | 13. Dieldrin |
| 5. / | 14. O,p'-DDT |
| 6. Methyl Parathion | 15. Beta-Endosulfan |
| 7. Malathion | 16. P,p'-DDT |
| 8. Fenotrothion | 17. / |
| 9. Ethyl parathion | 18. Tetradifon |

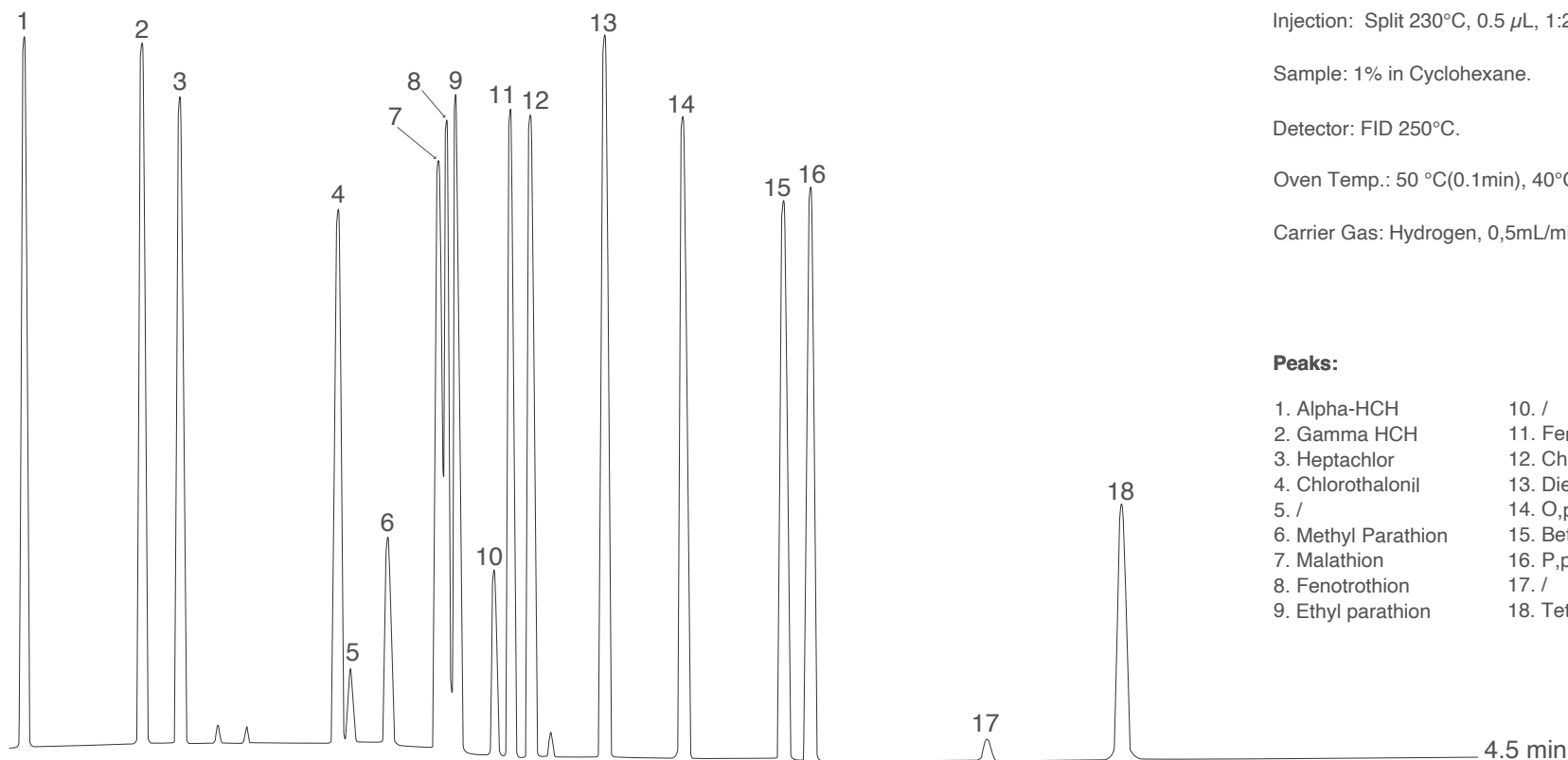

Acknowledgement: Prof C. Bicchi, C. Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.