


DL-Lactide

Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DEX DET Beta - 0.25mm, 0.25 μ m, 15m**
Catalog Code: DEX-DET-B-025-025-15

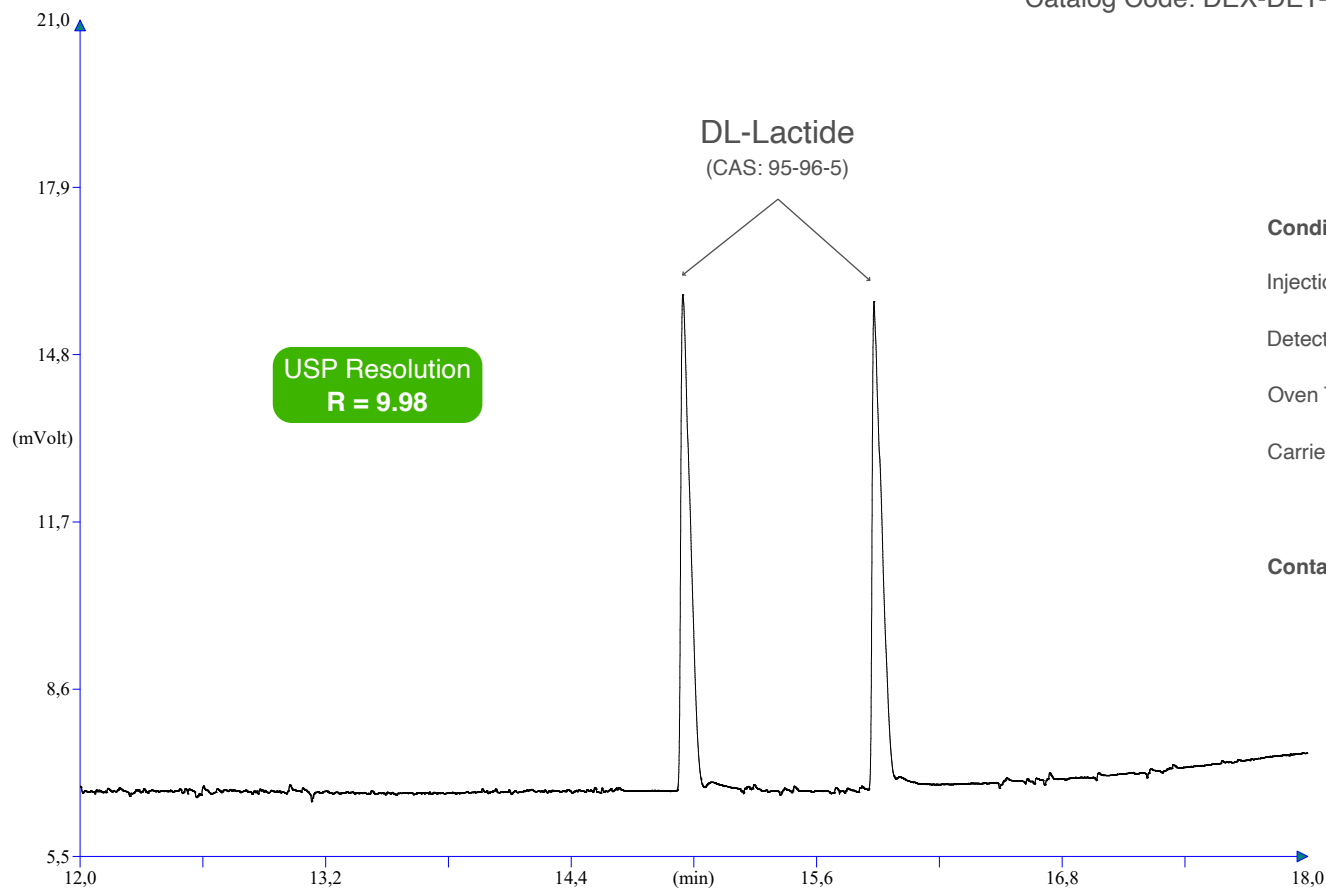

Conditions:

Injection: Split 1:100, 250°C, 0.5 μ L injected volume.

Detector: FID, 250°C.

Oven Temp.: 50°C, 5°C/min.

Carrier Gas: Hydrogen, 50 kPa.

Contact us for additional information.