

Ibuprofen (+/-)


Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA[®]
gas chromatography solutions

Column: **MEGA-DMP Beta 1701 - 0.25mm, 0.25 μ m, 15m**

Catalog Code: DEX-DMP-B-1701-025-025-15

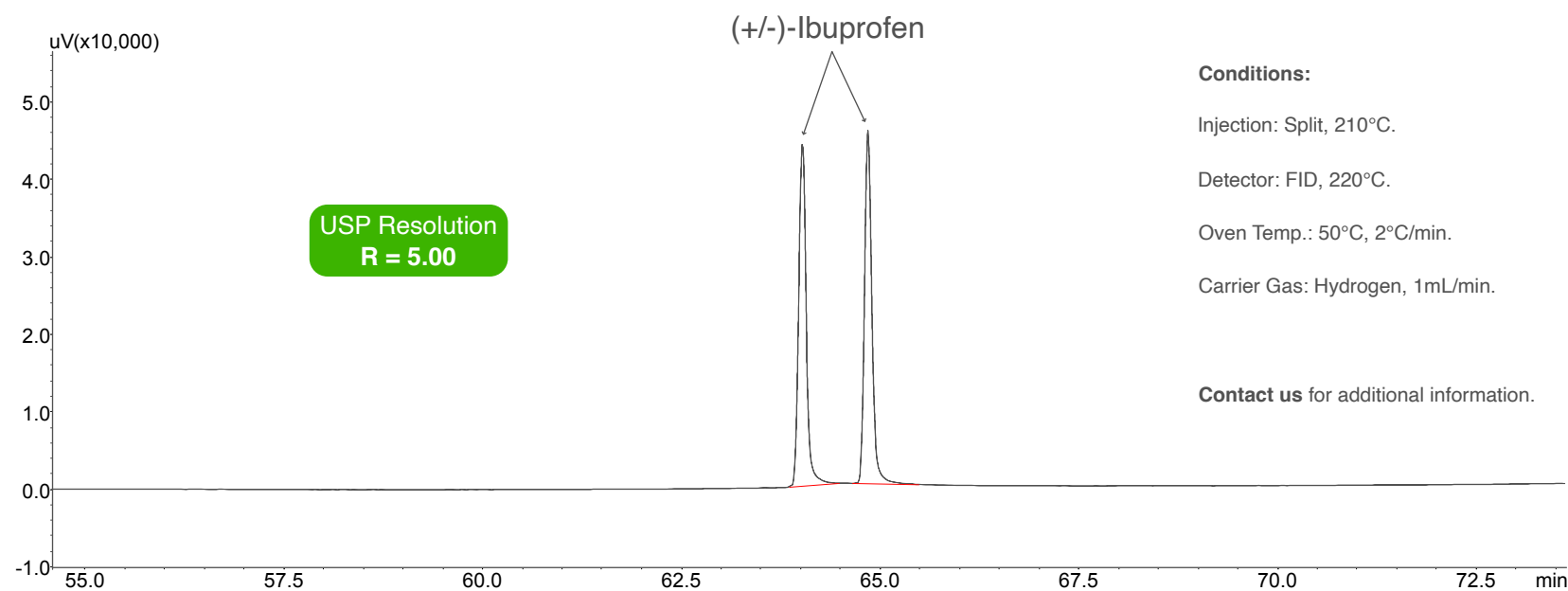

Conditions:

Injection: Split, 210°C.

Detector: FID, 220°C.

Oven Temp.: 50°C, 2°C/min.

Carrier Gas: Hydrogen, 1mL/min.

Contact us for additional information.

Acknowledgments: many thanks to Cecilia Cagliero, Gaia Bechis and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.