

## *cis/trans* - Nerolidol (+/-)

Chiral separation on MEGA-DEX columns

CHIRAL DEX

MEGA<sup>®</sup>  
gas chromatography solutions


Column: **MEGA-DEX DET Beta - 0.25mm, 0.25 $\mu$ m, 25m**

Catalog Code: DEX-DET-B-025-025-25

### Conditions:

Oven program: 50°C, 2°C/min.

Carrier Gas: Helium, 1mL/min.

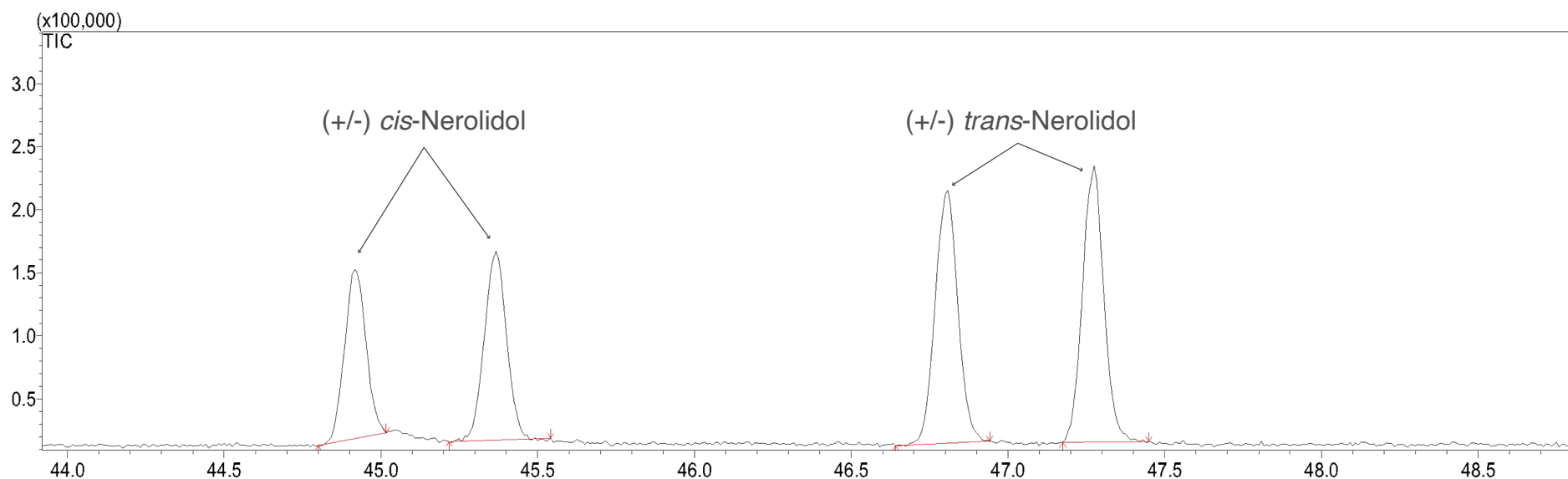

Acknowledgments: many thanks to Barbara Sgorbini, Cecilia Cagliari, Chiara Cordero and Carlo Bicchi, Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via P. Giuria, 9, Torino, Italy.